



**Powerology**

**Ultra-Mini Travel Massager**



SKU: PSM011

# Table of Contents

Specifications	2
Schematic View	3
Product Accessories	4
Instructions for Use	4
Operating Gears	5
Charging Instruction	6
Safety Precautions	6
Usage Prohibitions	7
Disposal	7
Warranty	8
Contact Us	8


Specifications

Model Number	SM011
Battery Capacity	750mAh
Recharge Time	1.5H
Working Time	Up to 4.5H
Maximum Power	40W
Charging Interface	Type-C
Gear Speed	1600 to 2800RPM
Stroke	4.5mm
Product Weight	200g
Product Dimension	93.5 x 76.3×35mm
Noise Level	About 50dB
Maximum Current	5A

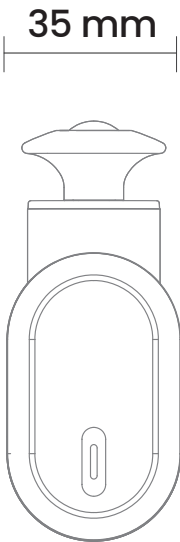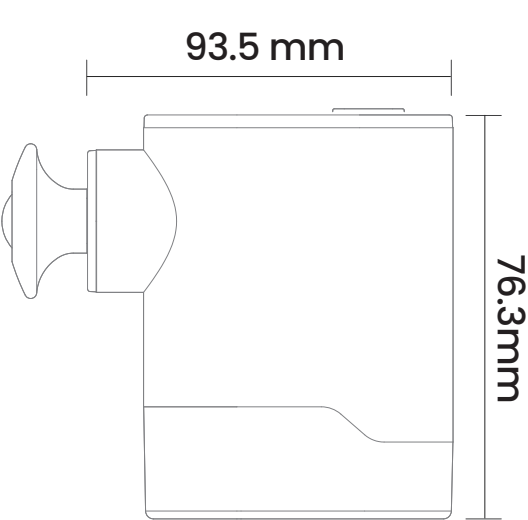

# Schematic View

1. Quick-Release Metal Massage Head
2. Grip
3. Massage Head Rubber Cover
4. Gear Indicator
5. Switch/Shift Compound Key/Power Indicator
6. Type-C Charging Interface

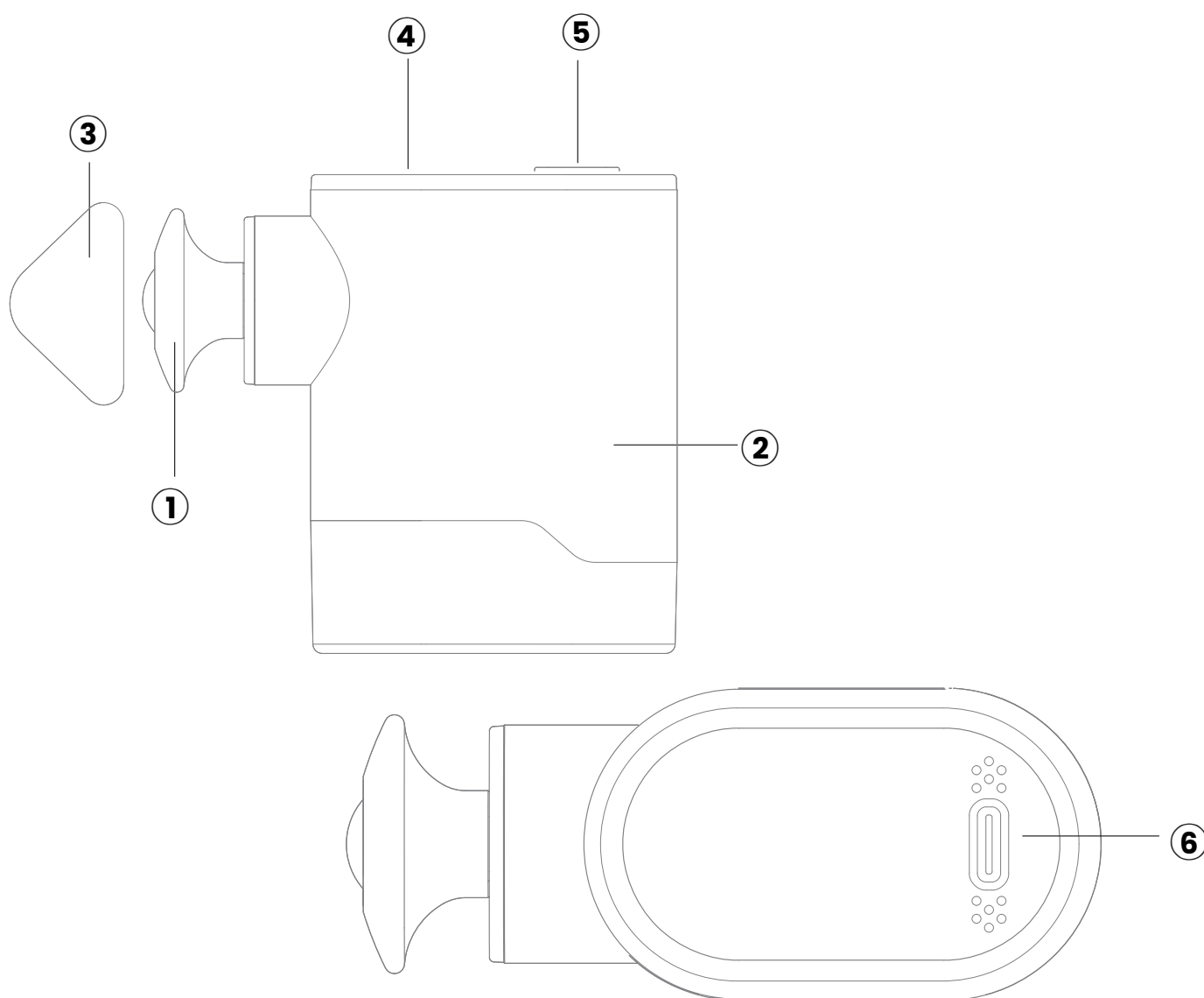

## Product Accessories

1. Flat Metal Massage Head
2. Triangular Massage Head
3. Hamburger-Shaped Massage Head
4. Straw Hat-Shaped Massage Head
5. Type-C Charging Cable



①




②



③



④

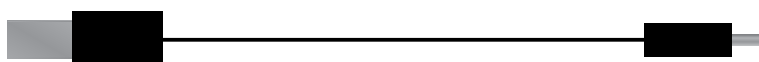

⑤

## Instructions for Use

### 1. Massage Head Installation

1.1 Please ensure the machine is turned off before installing the massage head.

1.2 Choose the appropriate massage head.

1.3 Gently insert the massage head in place.

## **2. Turning On/Off**

**2.1** To turn on: Long press the switch for about 1.2 seconds to enter gear 1; the green LED indicator will light up.

**2.2** To turn off: While the device is on, long press the switch for about 1.2 seconds until all indicator lights turn off, indicating that the device has entered shutdown mode.

## **3. Blocking Function**

**3.1** If the current exceeds 5 A, the device will stop working and enter standby mode with the gear indicator light remaining active.

**3.2** To resume operation from the current gear, short press the switch; the device will restart about 1 minute after an automatic shutdown.

## **4. Timing Function**

**4.1** The device runs for about 10 minutes and then automatically shuts down using a timer.

**4.2** To stop the automatic shutdown, press the switch once for 1 minute and again for 10 minutes.

## **5. Shift Control Function**

**5.1** To start in gear 1, short press the switch to increase through the gears.

**5.2** Press the switch three times to cycle back to gear 1.

**5.3** Gear indicators: 1st gear (Green), 2nd gear (Red), 3rd gear (White).

# **Operating Gears**

The machine starts in an initial state. To switch between modes, touch the shift switch to cycle through the gears from low to high, and neutral.

**Gear 1600** :1 rpm, white LED light – Suitable for warming up muscles.

**Gear 2200** :2 rpm, red LED light – Suitable for fascia relaxation.

**Gear 2800** :3 rpm, orange LED light – Suitable for lactic acid removal.

## **Charging Instruction**

1. When the power is insufficient, the power indicator will flash red 5 times and then turn off.
2. While charging, the power indicator light will remain red.
3. After charging is complete, the power indicator light will turn blue.

## **Safety Precautions**

1. This device is compatible with 5V / 2A chargers or partial fast charging adapters.
2. Fully charge the battery before using the device for the first time.
3. Ensure that the input voltage matches the adapter voltage before charging.
4. Do not change the massage attachments while the device is operating.
5. Do not immerse the product in water.
6. Avoid applying vibration to the head, groin, chest, or any areas where surgery has been performed.
7. Do not use the device as a substitute for medical treatment or therapy. If you have any pain or are undergoing treatment or using

other medical devices, consult a specialist before using this device.

**8.** In case of malfunctions, contact a technician. Do not attempt to disassemble the product yourself.

## Usage Prohibitions

Consult an expert before using the product if you have any of the following conditions or health problems:

- 1.** Pregnancy, diabetes, and complications such as nerve or retinal damage.
- 2.** The presence of a pacemaker, recent surgery, epilepsy or migraine, lumbar intervertebral disc prolapse, or vertebral displacement.
- 3.** Recently undergone joint replacement surgery, IUD placement, metal implants, or have other health concerns.

## Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.





## **Warranty**

Products that you buy directly from our **Powerology** website or shop come with a 24-month warranty.

The 24-month warranty applies to products purchased directly from our **Powerology** website or store. If **Powerology** products are bought from any of our verified retailers, then the product is eligible for only a 12-month warranty. To extend your product's warranty, visit our website **powerology.me/warranty** and fill in your details in the provided form along with an uploaded picture of the product to process your request. Once approved, you will receive a confirmation email of the extended product warranty. Upload the required information within 48 hours of purchase to be eligible for a 24-month warranty period.

For more info, please check:  
**powerology.me/warranty**

## **Contact Us**

If you have any questions about this Privacy Policy, please contact us at: **hey@powerology.me**

Website: **powerology.me**

Instagram: **powerology\_official**

Facebook: **powerology.ME**