

porodo



Porodo LifeStyle Lumbar Massage Pillow

SKU: PD-LSLMSP-GY

Table of Contents

Features	2
Specifications	2
Safety Guidelines	3
Usage Prohibitions	3
Warnings	4
Operating Environment	5
Product Disposal Safety	5
Battery Disassembly	5
Schematic View	6
Button Instructions	7
Troubleshooting	8
Environmental Notes	8
Maintenance and Storage	9
Warranty	10
Contact Us	10

Before installing and using the product, please carefully read this User Manual to guarantee correct usage and keep it secure for future reference.

Features

- 1. Headrest Area: Ideal for supporting the head during table naps.
- 2. Big Massage Nodes: Designed to target broader muscle regions.
- 3. Small Massage Heads: Intended for concentrated, deep pressure.
- 4. Shiatsu Kneading: Utilizes deep, circular motions and varying degrees of pressure to relax and ease muscle tension.
- 5. Hot Compress: Reduces stiffness in the muscles and provides a soothing experience.
- 6. Ergonomic Chair Strap: Easily secures to any chair for targeted lower back support.

Specifications

Battery Capacity	2600mAh
Rated Charging Voltage	DC 5V
Rated Power Input	10W
Auto Timer	15 minutes
Recharge Time	3-4H
Total Working Time	150 minutes (10 Minutes/session)
Product Net Weight	1.35KGS
Executive Standards	GB4706.1-2005, GB4706.10-2008

Note: This product is rechargeable and is compatible with Type-C interface charging

Safety Guidelines

Please read all the instructions and warnings in this manual carefully before using the product, to reduce the risks of failure, electric shock, injury, fire, death, and equipment damage.

1. If the machine does not work completely, stop using it immediately and turn off the power supply.
2. Do not operate the control panel with wet hands.
3. If using the product makes you feel uncomfortable, please discontinue its use.
4. If a power outage occurs during charging, turn off the power.
5. Please turn off the power when the product is not in use.
6. If the power cable is damaged, it must be replaced by professionals from the manufacturer's maintenance department or a similar department to prevent danger.
7. Do not use the product as a toy to prevent accidents.
8. Do not attempt to disassemble, modify, or fix the product yourself, nor unzip the inner sleeve yourself.
9. It is strictly prohibited to immerse the entire product in water for cleaning. If the product accidentally falls into water, it must be sent to a designated repair site; otherwise, it should not be used again.
10. Please keep the product away from sharp objects. If the outer cover is damaged, the product must not be used and should be replaced or repaired at authorized service centers.

Usage Prohibitions

If you have any of the following conditions, please consult your doctor before using the product:

1. Individuals currently receiving medical treatment or feeling unwell.
2. Patients diagnosed with diabetes, hypertension, heart disease, cerebral arteriosclerosis, or cerebral thrombosis.
3. Individuals diagnosed with malignant tumors.

4. Pregnant women or women who are in the process of giving birth.
5. Individuals with a body temperature exceeding 38°C.
6. Individuals insensitive to changes in skin temperature.
7. Individuals diagnosed with osteoporosis.
8. Individuals using pacemakers or equipped with medical implants.
9. Individuals with acute medical conditions.
10. Minors and individuals with physical and mental disabilities.

Warnings

1. This product is not a toy. Ensure children do not play with this product.
2. This product is not recommended for use by individuals (including children) with physical disabilities, sensory or mental disorders, or lack of relevant experience and knowledge, unless supervised and guided by responsible personnel.
3. As this product includes heating devices and the surface can become very hot, individuals who are not sensitive to heat should exercise caution when using it.
4. If the power cord is damaged, it must be replaced by the manufacturer or a designated maintenance unit or professional maintenance personnel to avoid hazards.
5. Do not apply excessive force or use fingernails on the operation panel. Please gently touch the switch.
6. Do not let the control panel come into contact with liquids such as water, and it is forbidden to operate the control switch with wet hands.
7. Do not insert or drop any items into the openings of the product.
8. Do not fall asleep while using the product.
9. Avoid using the product within one hour before or after eating to prevent adverse effects.
10. Do not use this product if you are intoxicated or feeling unwell.
11. If you feel uncomfortable during use, discontinue use immediately.
12. This product is timed for 10-minute sessions. Do not use it for more than 10 minutes at a time, and allow the device to cool down for 30 minutes after use to extend its service life.
13. Please turn off the power switch of the product in a timely manner after use.
14. If the device malfunctions, STOP using it immediately.

15. This product should not be used for self-diagnosis or treatment. If necessary, consult a medical professional.

Operating Environment


- 1.** Do not use this product near the stove, inside electric heating quilts, or on heating appliances such as electric blankets, as it may cause a fire or malfunction.
- 2.** Do not cover the product with towels, cloths, or other objects to avoid overheating the internal parts and causing damage or failure.
- 3.** Avoid using the product in humid environments such as bathrooms or near water sources, as this may lead to malfunctions and endanger personal safety.
- 4.** Do not use the product outdoors or expose it to direct sunlight for extended periods.

Product Disposal Safety


This product should not be discarded with household waste. Before disposal, please follow these steps to remove the lithium battery and give it to a professional organization capable of recycling batteries, or place it into the appropriate recycling equipment.

Battery Disassembly


- 1.** Before removing the battery, ensure that the product is powered off.
- 2.** Remove the fabric cover of the product and extract the movement circuit board.
- 3.** Remove the cover of the movement circuit board box.
- 4.** Disconnect the battery by unplugging the battery cable plug.
- 5.** Hand over the dismantled battery to a professional recycling agency and use professional recycling equipment.

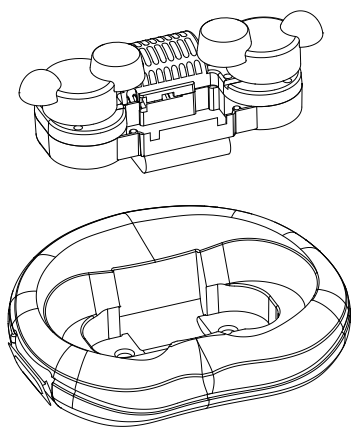

Step 1

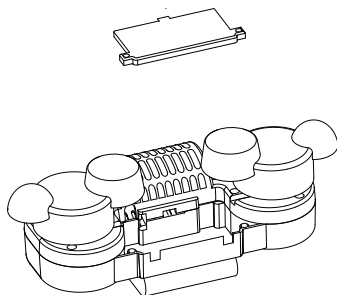


Step 2

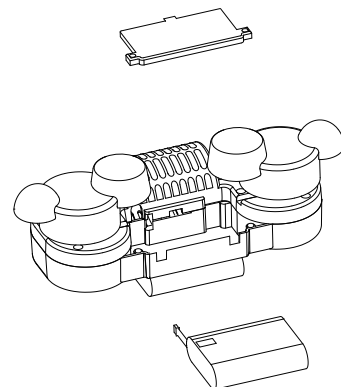


Step 3


Schematic View

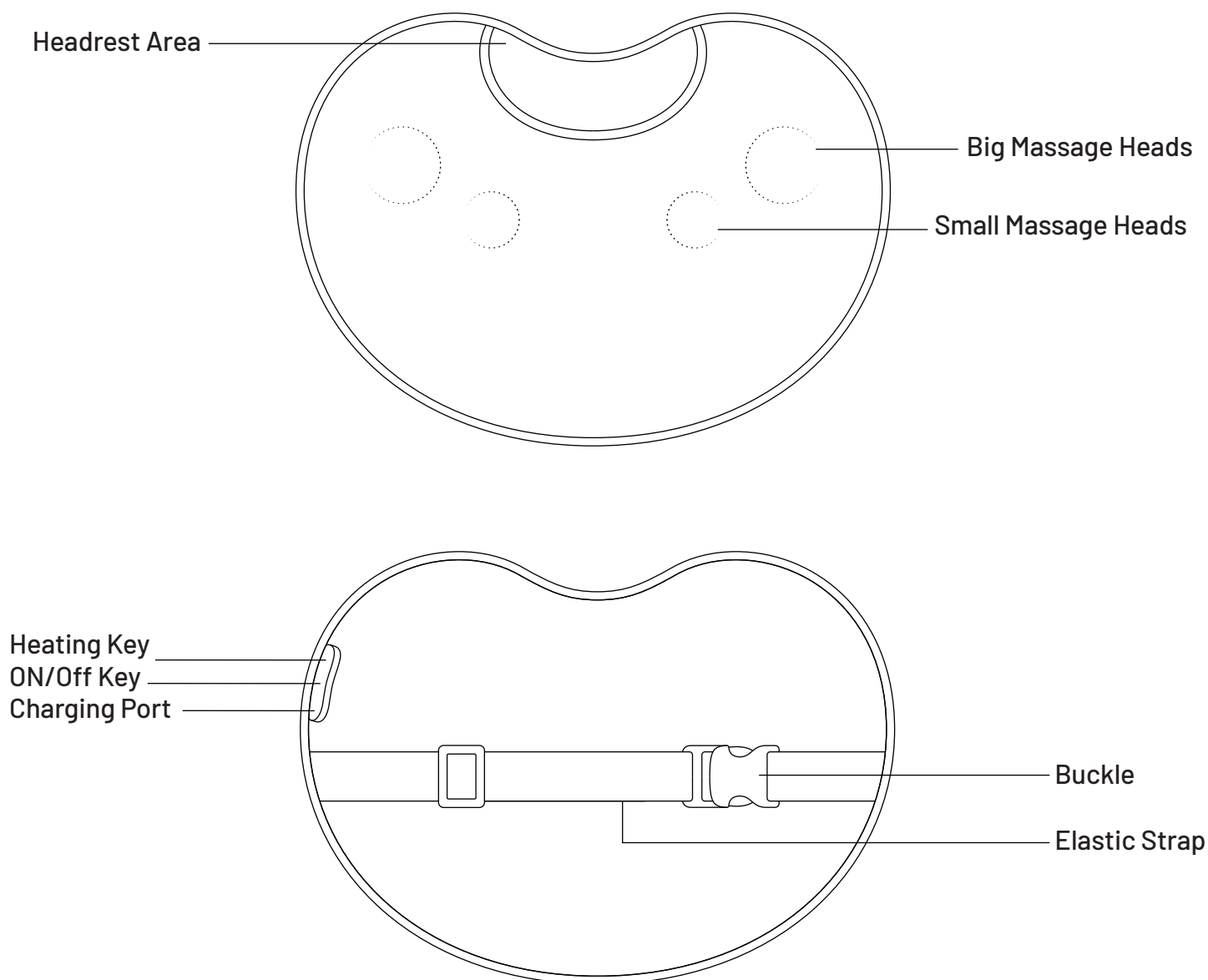


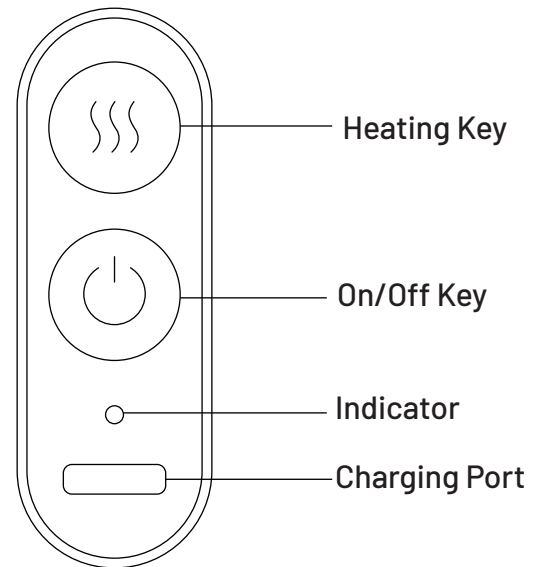
Button Instructions

On/Off Key


1. Press and hold the  button for 2 seconds; the buzzer will sound, and the device will begin operating. The default operating time is 15 minutes.

2. To turn off the device, press and hold the  button again for 2 seconds. The buzzer will sound, indicating that the device has stopped working.

3. In the active state of the device (default mode one), briefly press the  button; the buzzer will sound once, allowing you to switch between mode two and mode three and cycle through the settings accordingly.





Heating Key



While the device is operating, briefly press the  button to toggle the hot compress function on or off. A buzzer will sound to indicate the change.

Note: Pressing the  button in mode three is not functional.

Low Battery Reminder

When the battery level is low, the buzzer will sound twice to signal a low battery. Simultaneously, the indicator light will flash red quickly  . Please ensure the device is charged promptly.

Charging Indicator

The indicator light remains red during charging . When the device is fully charged and in operation, the indicator light will emit a solid green light .

Troubleshooting

Fault Phenomenon	Possible Cause	Solution
Unable to start	1. The device is not powered on. 2. Check if the battery voltage is low. 3. Failure of internal parts.	1. Press and hold the On/Off Key. 2. Charge the device. 3. Contact the designated service center for the product.
Short battery life; Battery cannot be charged	1. The battery is degraded. 2. Failure of internal parts.	1. Send the device for battery replacement. 2. Send it to the maintenance site for inspection.
Abrupt stop during use	1. The battery level is low. 2. The timer has ended.	1. Charge the product in a timely manner and ensure it is fully charged before use. 2. Restart the product and continue using it.
Smoke and odor during operation	Internal part failure.	Immediately stop using the product and contact the Service Center.

Note: In the event of product malfunction or failure to start, please consult the aforementioned tips and troubleshooting methods. Should the issue persist, refrain from further use to prevent potential accidents and promptly reach out to authorized service and repair centers for assistance. Avoid disassembling or attempting to repair the product independently.

Environmental Notes

Please refrain from disposing of this product with regular household refuse once it reaches the end of its service life. Instead, it should be taken to a recycling facility sanctioned by the government. Staff at the disposal center will ensure that the product undergoes environmentally responsible processing.

To ensure commitment to environmental stewardship, follow these protocols:

1. Adhere to environmental laws and regulations.
2. Implement initiatives for conserving resources, optimizing energy use, and minimizing waste across all company operations.

3. Strive for the reduction of environmental pollution and foster the enhancement of our local surroundings.
4. Integrate environmental protection efforts with our business growth strategies.
5. Motivate our workforce to actively engage in measures that positively impact the environment.

When disposing of our product packaging, please adhere to the following guidelines:

1. The outer packaging carton should be processed according to local recycling regulations.
2. Expandable Polystyrene, characterized as solid form, is recyclable. You can place this material in your local recycling bin.
3. Plastic packing bags are also recyclable and should be deposited in the appropriate recycling container.
4. For packing bags, follow local guidelines to determine the proper disposal method.
5. Lastly, bubble packs require specific handling; consult with local waste management services to ensure correct disposal or recycling.

Maintenance and Storage

A. Storage

1. If the product is not in use for a long period, please fully charge it and store it in a safe, dry, and cool place.
2. It is recommended to fully charge the battery every two months to prevent battery damage due to power depletion.
3. Avoid using sharp or hard objects that could scratch the surface of the product.
4. Do not expose the product to direct sunlight or high temperatures.
5. Keep the product out of reach of children.

B. Maintenance

1. Ensure the power switch is turned off before cleaning.
2. Wipe the product with a soft and dry cloth.
3. If the product is particularly dirty, use a damp cloth with a neutral detergent, and gently clean the product.

4. Allow the product to dry completely after cleaning.

Warning: Do not use volatile oils, diluents, alcohol, strong alkaline detergents, or bleaching agents to clean the product.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **<https://porodo.net/warranty>** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:

<https://www.porodo.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **[porodo](https://www.instagram.com/porodo)**