



Green Lion

Mini Massage Gun 6 Speed Levels



SKU: GNMINMSGN

Model No: GL-MG28

Table of Contents

| | |
|----------------|---|
| Quick View | 2 |
| Usage Method | 2 |
| Warning | 4 |
| Specifications | 5 |
| Disposal | 6 |
| Warranty | 6 |
| Contact Us | 6 |

Quick View

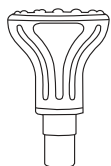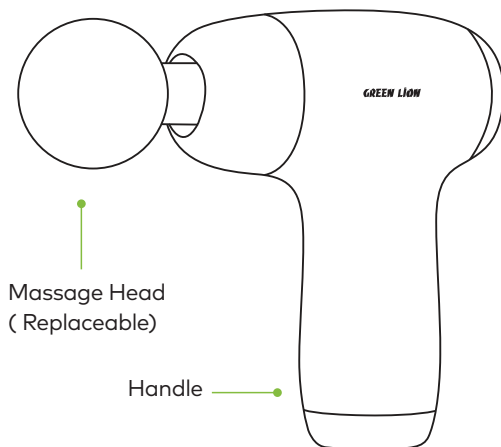


Flat Head



U-Shaped Head



Cylindrical Head

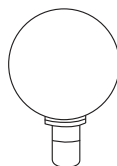

Ball Head

Usage Method


On/Off

1. Press the key to start the machine.
2. Press the key for two seconds to shut down.

Strength Adjustment

1. The initial gear on startup is the first gear.
2. Click the key once to increase by one gear.
3. The maximum gear is the sixth gear. When at the sixth gear, press and hold for 2 seconds to shut down.


Replacement of Massage Head

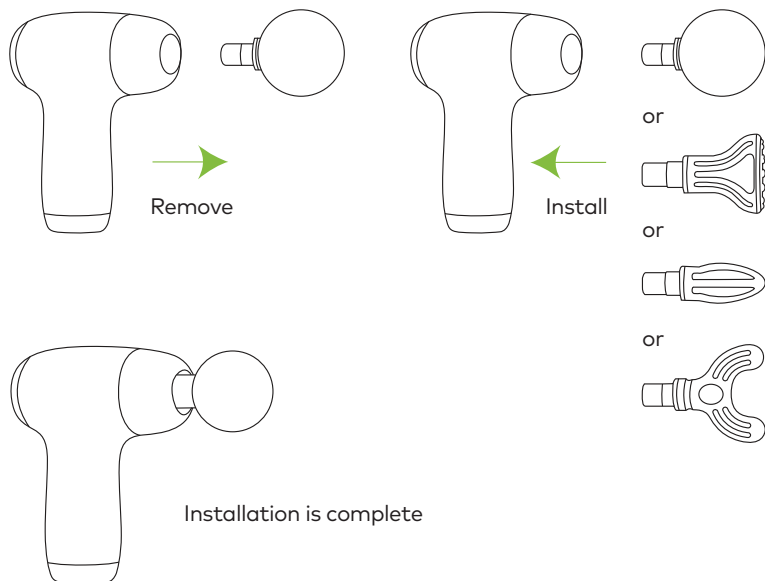

Indicator light

1. The indicator light shows the gear size of the fascia gun.
2. It is the first gear when one light is on.
3. The second gear is indicated when two lights are on.
4. The third gear is indicated when three lights are on.
5. The fourth gear is indicated when four lights are on.
6. The fifth gear is indicated when five lights are on.
7. The sixth gear is indicated when six lights are on.

Charge

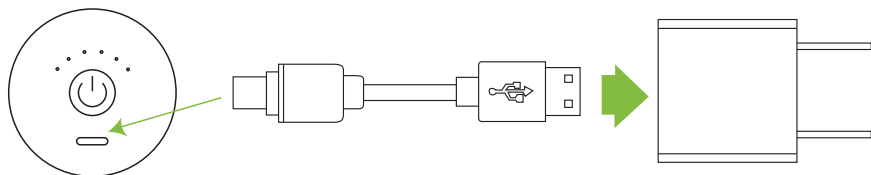

1. Please ensure that the battery is fully charged before using the device for the first time.
2. The indicator light flashes while charging and remains on when fully charged.
3. When the machine is off, press the power key briefly, and the indicator light will show the power status.

Indication

1. Relieve pain and spasms caused by muscle contusion, sprain, and strain.
2. Help the flow of edema fluid in the swollen area.
3. Loosen and reduce thickened connective tissue and fascia.
4. Decrease the accumulation of lactic acid in muscles.
5. Increase joint range of motion.
6. Eliminate muscle fatigue.

Forbidden

1. Aneurysms, bleeding, or use of blood thinners.
2. Heart disease, pacemaker, or defibrillator.
3. Pregnancy or cancer.
4. No 3-inch internal fixation is allowed within 90 days.
5. Use within the scope:
6. Sensitive areas such as the head, face, cervical spine, and vertebral body.
7. Muscle coverage that is less close to the bone, such as the tibia, dorsum of the foot, dorsum of the hand, and similar areas.

Warning

1. Do not drop or collide with the charger (adapter). If the shell is damaged, please contact the authorized service center for replacement.
2. Do not touch the power cord with wet hands or pull the power cord. Always pull out the charger (adapter).
3. Please do not touch the product or charger (adapter) with wet hands to avoid any damage, such as product short circuit, faults, or electric shocks.
4. Please stop using the charger if it is drenched by rain, liquid, or seriously affected by damp. Contact the authorized service center for replacement.

5. Do not expose the battery to high temperatures or around heating equipment, such as sunshine, heater, microwave, or water heater, etc., as overheating may cause explosion.
6. Do not disassemble or modify the battery, insert foreign materials, or immerse it in water or other liquids to avoid battery leakage, overheating, fire, or damage/explosion.
7. Before discarding the appliance, the battery must be removed from the appliance.
8. When removing the battery, the appliance must be powered off.
9. The battery should be disposed of safely.

Specifications

| | |
|--------------------|---------------|
| Model No. | GL-MG28 |
| Material | ABS |
| Motor | Brushed Motor |
| Battery | 1200mAh |
| Power | 24W |
| Voltage | 7.4V |
| Speed | 6 Levels |
| RPM Range | 1200–3200 |
| Charging Port | Type-C |
| Charging Time | < 2 hrs |
| Working Time | 1–2 hrs |
| Auto Shut-Off | 15 min |
| Product Weight | 500 g |
| Product Dimensions | 105 × 130 mm |

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Green Lion** website or shop come with a 24-month warranty.

When you buy **Green Lion** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at <https://www.greenlion.net/warranty> and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
<https://www.greenlion.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

Website: <https://www.greenlion.net/>

Sales: sales@greenlion.net

Service Support: Support@greenlion.net

Telegram: [@greenlion_global](https://t.me/greenlion_global)

Instagram: [@greenlion_global](https://www.instagram.com/greenlion_global)