

porodo



Porodo LifeStyle

Mini Massage Gun with 4 Head

SKU: PD-LSMSGN

Table of Contents

Main Objectives	2
Warnings	2
Schematic View	3
Specifications	3
Massage Head Functions	4
Massage Head Replacement	4
Operation Guide	4
Charging Guidelines	5
FAQ	5
Maintenance Guidelines	6
Environmental Protection	6
Warranty	7
Contact Us	7

To mitigate the risk of malfunction, electric shock, injury, fire, death, and equipment damage, it is crucial to thoroughly read all instructions and warnings provided in this manual.

Main Objectives

This device is primarily designed to alleviate muscle soreness, enhance blood circulation, and alleviate cervical stiffness.



Warnings

Prohibitions


1. Minors are strictly prohibited from using this product as it is not a toy.
2. In case the product accidentally falls into water, please promptly arrange for inspection.
3. Avoid placing the product near a heat source or exposing it to fire.
4. Refrain from using the product around the head, neck, heart, and genitals.
5. Only apply the product to soft tissues of the body; avoid using on bones.
6. Use only original chargers for charging to prevent potential fire hazards.
7. Intentional dropping of the product is prohibited.
8. Do not remove screws or attempt to disassemble the product.
9. Use of equipment as a substitute for medical treatment is prohibited; always consult a physician before using the device.
10. Keep fingers, hair, or any body parts away from the shaft and back of the massage head to prevent pinching.

When starting the machine without the massage head installed, avoid touching this area.

When installing the massage head and initiating the machine, refrain from touching this area.



Schematic View

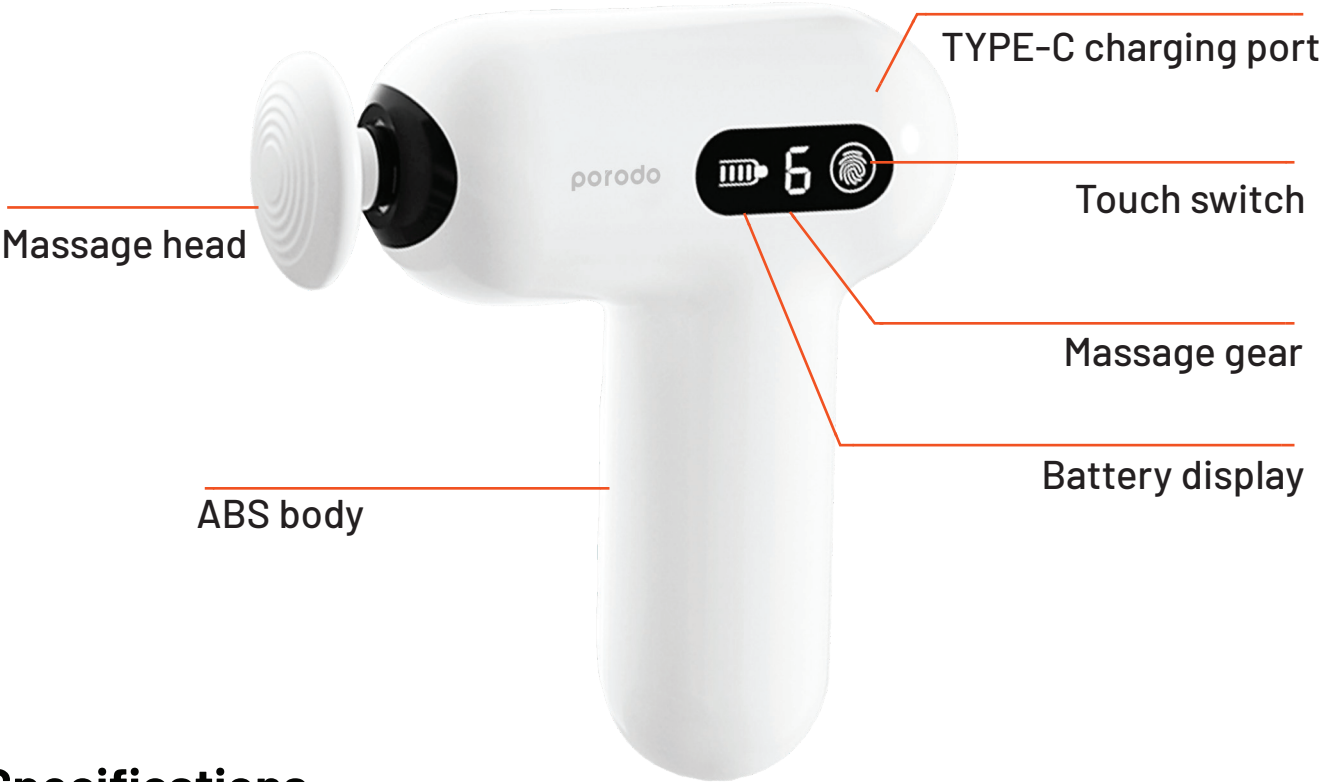

Specifications

Product Name	Mini Massage Gun with 4 Head
SKU	PD-LSMSGN
Interface Type	Type-C
Charging Input	5V=2A
Massage Stroke	6mm
Battery Capacity	1500mAh
Number of Gears	6
Button Method	Touch the button
Unit Weight	~ 433g (with T-shaped head)
Size	54×117×149mm
Operating Environment	0°C-45°C
Executive Standard	GB4706.1-2005 & GB4706.10-2008


Massage Head Functions



U-shaped head: Designed for spinal massage, providing intense stimulation.	Hemispherical head: Primarily utilized for muscles in the legs and buttocks.	T-shaped head: Ideal for relaxing large muscle groups, offering strong stimulation.	Bullet head: Targets deep joints for thorough relaxation.
---	---	---	--

Massage Head Replacement

1. Ensure the machine is powered off.
2. Carefully insert the massage head into the barrel of the massage gun.
3. Align the red line position of the massage head with the corresponding red line position on the machine, as illustrated in the figure.

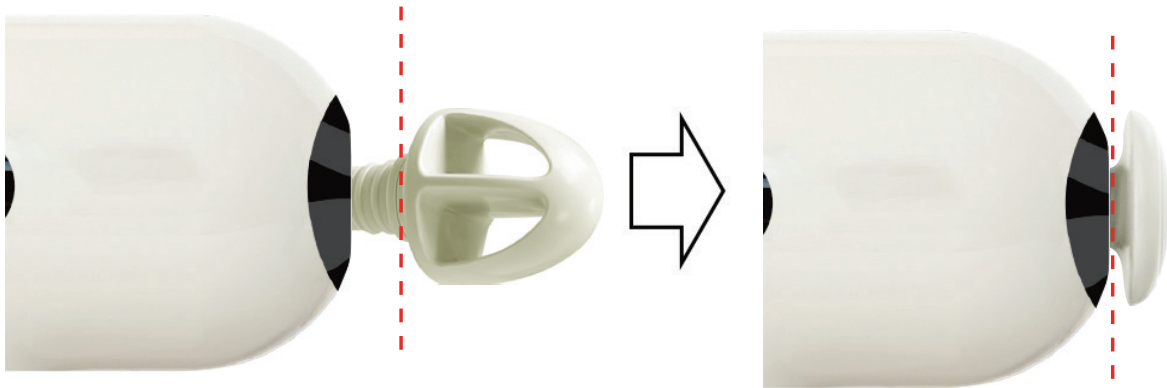

Operation Guide

1. Startup and Gear Switching: Press and hold the touch switch for more than one second to activate the device. Upon startup, the display will show 0 gear. Press briefly to switch to the first gear. Each subsequent short press increases the vibration intensity by one level, up to a maximum of 6 gears. The display screen will indicate the current gear number.

2. Shutdown: To turn off the device, press and hold the touch switch while it is in operation for one second.
3. Battery Indicator: The display screen indicates the current battery level. If the battery level is low, the battery icon will flash five times and then the device will shut down automatically to prompt charging.
4. Timing Function: The device will automatically stop after 10 minutes of continuous operation.
5. Locked Rotor Protection: If the device is obstructed for more than two seconds, it will enter a protection mode. The display screen will flash as an alert. Pressing the touch switch briefly within five seconds will restart the device. If there is no interaction, the device will automatically power down.

Charging Guidelines

1. Before initial use, ensure the machine battery is fully charged.
2. Switch off the machine and connect the Type-C charging cable to the charging head (ideally 5V2A). The battery icon will flash, indicating charging initiation.
3. Once the battery icon remains steady, it signifies a complete charge.
4. Avoid operating the machine while it's charging; it cannot be turned on during this process.

FAQ

Issue: The device won't Turn On

Possible Cause 1: Low Temperature

Extended exposure of the fascia gun to temperatures below 0 degrees Celsius can diminish battery performance, causing startup issues. While the machine may appear to initiate normally, it fails to operate effectively in either direction.

Possible Cause 2: Insufficient Battery Power

If left uncharged for an extended period, the battery enters a low-power protection state, rendering the fascia gun inoperable. Solution: Ensure normal functionality by fully charging the battery.

Operational Troubleshooting	
Product Non-Functionality	Verify the adequacy of power supply.
Panel Light Inactivity	Confirm the sufficiency of power supply.
New Machine Emitting Sound	This is considered a normal occurrence.
Battery Charging Failure	If the battery is aging, replacement is necessary.

Maintenance Guidelines

- 1. Long-Term Storage:** When not in use for an extended period, maintain the battery at full capacity, switch off the power, and store it in a safe, dry, well-ventilated, and cool environment.
- 2. Periodic Charging:** Charge the product every six months when not in use for an extended duration.
- 3. Daily Care:** Utilize a neutral detergent and a soft cloth for regular maintenance.
- 4. Avoid Chemicals:** Refrain from using oil, organic solvents, or other chemical cleaning agents. Avoid rinsing with water.
- 5. Power Management:** Always switch off the power before using or cleaning the product.
- 6. Secure Transportation:** Ensure the product is securely fixed during transportation to prevent rolling or falling.

Environmental Protection

At the conclusion of the product's service life, dispose of it in Government-designated recycle bins to ensure environmentally responsible recycling.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **<https://porodo.net/warranty>** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:

<https://www.porodo.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **[porodo](https://www.instagram.com/porodo)**