

PORODO
LIFESTYLE



Porodo Lifestyle
All-Purpose Massager
with Red Light and Hot Compress

SKU: PD-LSNBSM-WH

Table of Contents

Features	2
Safety Warning	2
Battery Disposal	2
Usage Prohibitions	3
Charging Instructions	3
Instructions for Use	3
Warranty	5
Contact Us	5

Features

1. Offers two distinct massage settings to cater to varying relaxation needs and preferences.
2. Features a simulated kneading function that mimics the pressure and movement of a professional massage, providing deep muscle relief.
3. Equipped with a hot compress feature that applies gentle heat to soothe and relax tight muscles, enhancing the therapeutic effects of the massage.

Safety Warning

1. Before using the product, please carefully read this User Manual to guarantee correct usage and keep it secure for future reference.
2. If the appliance ceases to function or exhibits abnormal behavior, discontinue use immediately and seek professional maintenance. Do not attempt repairs yourself.
3. Limit use to no more than 30 minutes at a time to prevent overheating or discomfort.
4. When disconnecting the power, grasp the plug rather than pulling on the cord to avoid damage.
5. Do not use the appliance in bathrooms or areas with high humidity. Avoid placing it near water sources to prevent the risk of electric shock or short circuit.
6. This appliance is not suitable for children and should not be used by individuals who are not sensitive to heat due to its heating surface.
7. Remove the battery from the appliance before disposal of the product.

Battery Disposal

1. Ensure the appliance is turned off before removing the battery.
2. Authorized personnel must use a screwdriver to open the casing of the product. Once the casing is removed, the battery can then be safely extracted.


3. Dispose of the battery in accordance with local environmental regulations to ensure safe disposal.

Usage Prohibitions

1. Individuals undergoing treatment or with certain medical conditions should consult a physician before using this appliance. This includes those with heart disease, malignant tumors, those who are pregnant, individuals with back abnormalities, osteoporosis, or those with an implanted cardiac pacemaker.
2. This appliance is not suitable for individuals, including children, who have physical or sensory impairments, intellectual disabilities, or lack experience and knowledge concerning its operation. Such individuals should only use this appliance under the supervision and guidance of someone responsible for their safety. Ensure the appliance is turned off after use.

Charging Instructions

1. Ensure the battery is fully charged before using the device for the first time.
2. While charging, the indicator light will display red. Once the battery is fully charged, the light will change to green.

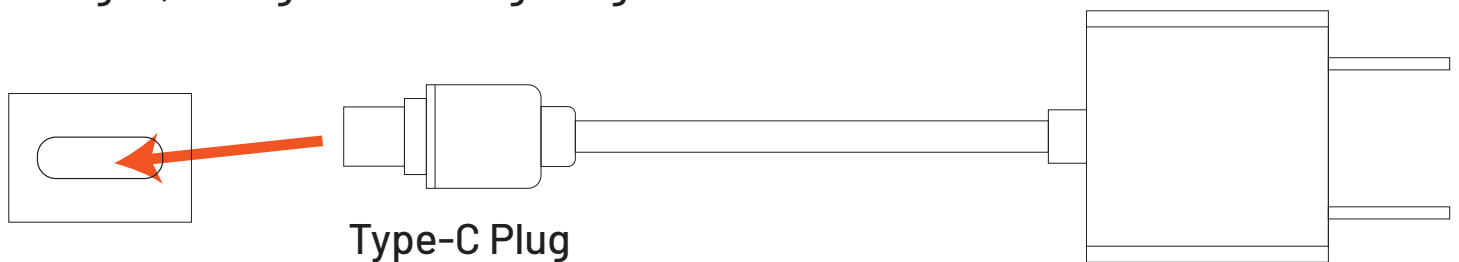

Instructions for Use

1. Before use, inspect the massager for any damage. To protect the skin, apply a generous amount of essential oil to the area to be massaged.
2. Place the massager on the area to be treated.
3. Press the "power switch" key on the control panel to start the device. By default, the massager runs for 10 minutes per session.
4. Press the "direction" key to alternate the rotation of the massage head between

- clockwise and counterclockwise to vary the massage effects.
5. Use the "speed regulation" key to change the operating speed of the massager according to your preference.
6. Activate or deactivate the light wave function by pressing the "light wave" key.



Specifications

Product Name	Porodo Lifestyle All-Purpose Massager with Red Light and Hot Compress
SKU	PD-LSNBSM-WH
Charging Current	2A
Rated Power	24W
Battery Capacity	2000mAh
Size	10×15.6×26.5cm
Net Weight	1.15kg

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**