

PORODO
LIFESTYLE



Porodo Lifestyle
4 In 1 RF, EMS, Vibration Light
Therapy Device Face

SKU: PD-LFST100

Table of Contents

Features	2
Specifications	2
Usage Prohibition	2
Warnings	3
Important Safety Instructions	4
Schematic View	5
Disposal	5
Warranty	6
Contact Us	6

Features

- 1. Strong RF Activation:** Stimulates collagen production, enhances elasticity, and improves blood circulation for a radiant complexion.
- 2. EMS Micro-Stimulation:** Stimulates muscles, promotes blood flow, and revitalizes the skin, promoting a toned and firm appearance.
- 3. Vibration Massage:** Relieves tension, boosts blood circulation, and leaves the skin feeling refreshed, relaxed, and visibly energized.
- 4. LED Light Therapy:** Reduces wrinkles, clears acne, and repairs skin using red blue, and green light wavelengths for optimal skin health.

Specifications

Battery Capacity	1450mAh
Power	6.5W
Voltage	5V/1A
Charging Interface	Type-C
Modes	RF Mode, EMS Mode, VR Mode, LED Mode
Vibration Frequency	1500±8500rpm
Product Size	180×44×44mm
Product Weight	118g

Usage Prohibition

This product should not be used by the following individuals (as it may cause accidents or injuries):

- 1. People with artificial cardiopulmonary devices or life-sustaining medical**

equipment.

2. People using medical equipment such as pacemakers.
3. Patients with visceral diseases or facial nerve disorders.
4. Individuals with mental illnesses or acute diseases.
5. People with abnormal blood pressure or infectious diseases.
6. Individuals with tuberculosis or blood disorders.
7. Pregnant women, menstruating women, breastfeeding mothers, or those with a fever.
8. Those undergoing sexual treatment or with atopic skin conditions.
9. People with allergies or those with heightened skin sensitivity.
10. Individuals consuming alcohol, taking sedative drugs, or experiencing excessive fatigue.
11. Long-term users of steroid hormones or patients with telangiectasia caused by liver dysfunction.

Warnings

The following areas should not be used (otherwise accidents or injuries may occur):

1. Wounds, plastic surgery sites, eyeballs, and upper eyelids.
2. Inside the mouth, pubic area, mucous membranes, etc.
3. Areas with excessive burns, skin itching, swelling, inflammation, pain, or perceptual disturbances.
4. Areas where metal, plastic, silicone, etc., are embedded under the skin.
5. In LED mode, the light can shine around the eyes. When using, close your eyes and wear an eye mask.
6. Do not look directly at the LED light.
7. Do not stay in the same area for too long as it can cause burns or frostbite.
8. Do not charge or use the device in a fire or high-temperature environment, as this can lead to explosion or rupture of the internal battery.
9. Keep out of reach of children, as this could lead to accidents or injuries.
10. Follow the instructions for proper use. Not doing so may result in injury or

machine failure.

11. If unsure, consult a doctor before use.

12. If you experience pain, discomfort, or abnormal skin reactions after use, stop immediately and consult a doctor.

13. First-time users should adjust based on their skin condition. Failure to do so may cause skin damage or discomfort.

14. Remove any rings, necklaces, or other accessories, and avoid wearing

Important Safety Instructions

1. Do not use the device while taking a shower or immerse it in water.

2. If any abnormality occurs during charging or use, immediately disconnect the power plug from the charging socket.

3. Only use power supplies with specifications of AC 240-100V and DC 5V 1A. Using other power supplies may damage the device.

4. Do not use power plugs that do not match the charging voltage required by the device. Using the wrong plugs can damage the device.


5. Do not use damaged power cords, plugs, or sockets. These could cause a fire or electric shock.

6. Do not charge the device in a shower or damp environment. This can result in electric shock or short-circuiting.

7. Do not insert or remove the charging plug with wet hands. Doing so can cause device failure or a short-circuit.



Schematic View

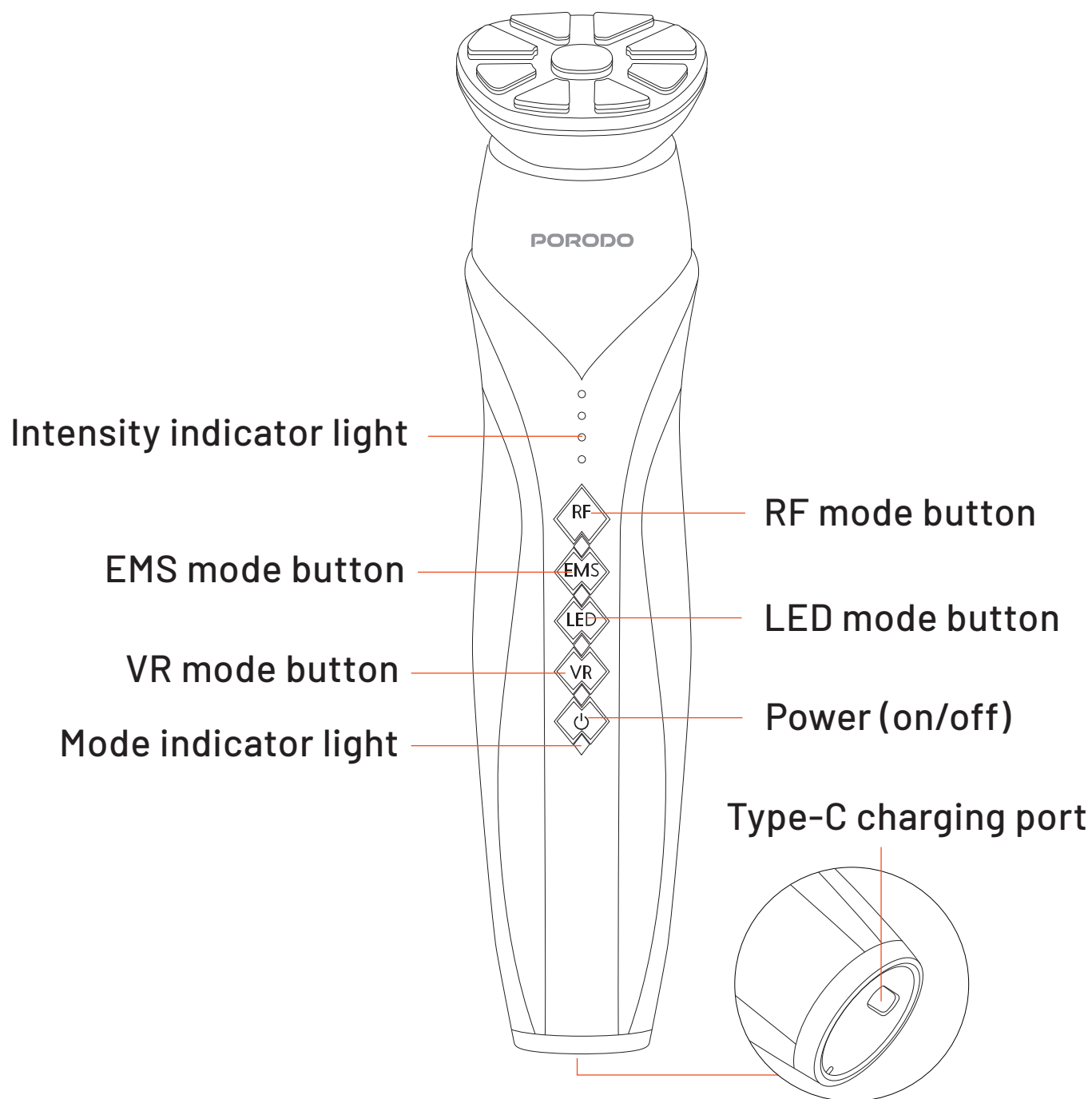

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**