

PORODO
LIFESTYLE



Porodo Lifestyle


IPL Hair Removal

SKU: PD-LFST113


Table of Contents

Package Contents	2
Schematic View	2
Product Overview	2
IPL Principle of Operation	3
Specifications	4
Instructions for Use	5
Applicable Skin and Hair Color	6
Applicable Body Parts	6
Troubleshooting	7
FAQ	9
Disposal	10
Warranty	11
Contact Us	11

Package Contents

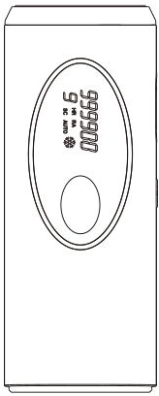


Hair Removal Unit

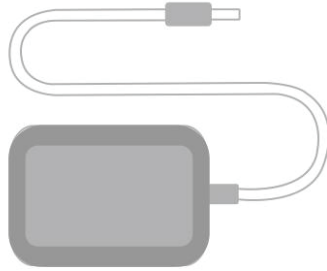

Adapter




Goggles

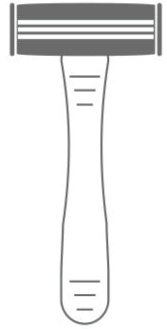

Razor


Schematic View

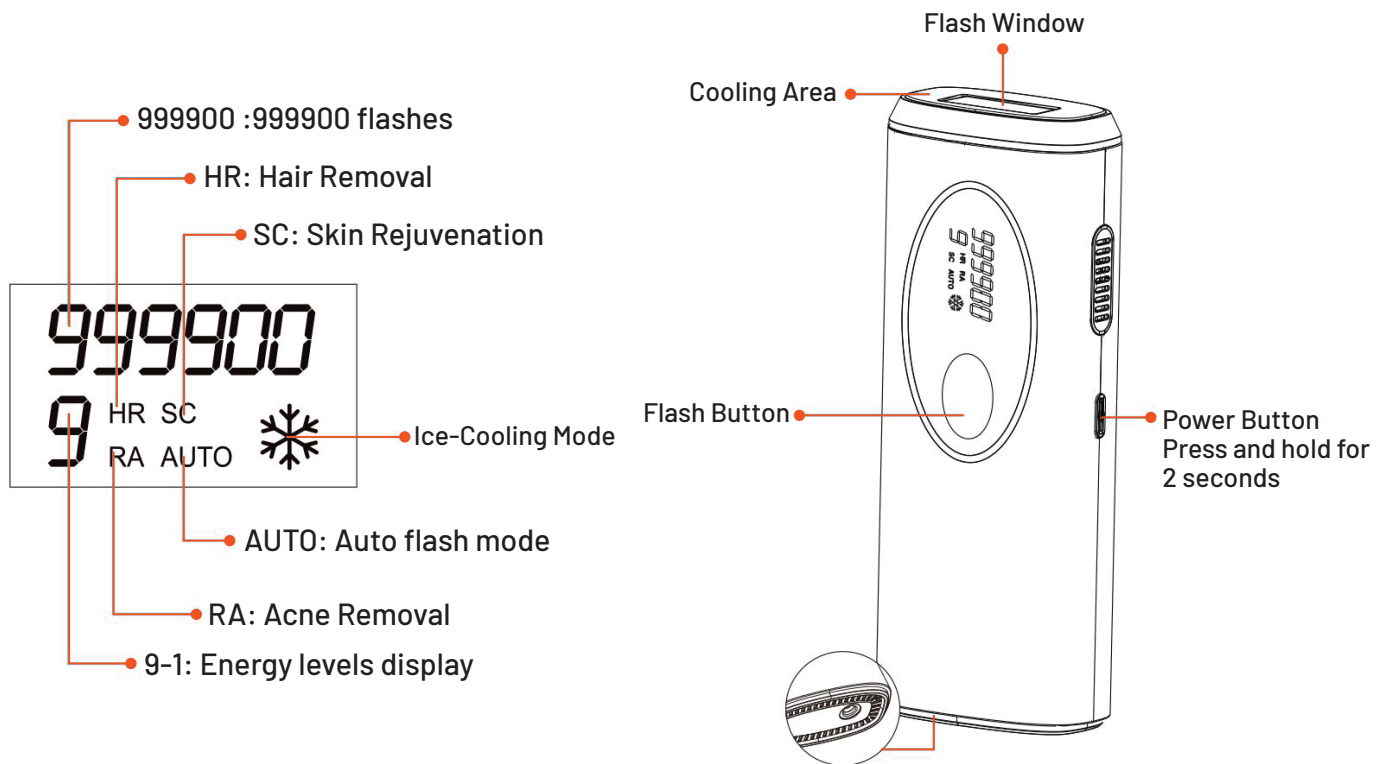

Product Overview


Under the exposure to pulsed light with a specific wavelength and pulse width the melanin in the skin absorbs energy and transfers the transient heat to the hair follicle. As the hair follicle reaches a certain temperature, it enters a dormant

phase, causing the hair to stop growing and eventually fall out. After using the device for about 8 weeks, the desired results can be achieved.

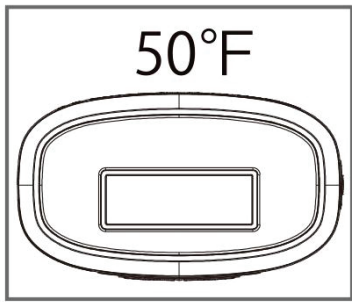

Enabling the ice-cooling mode helps to alleviate burning and irritation during operation, providing a painless experience.

Additionally, the flashing speed has been increased allowing you to free your hand for a more efficient treatment.

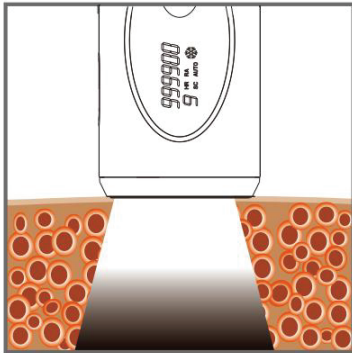


The multi-functional hair removal device is an innovative beauty product designed to effectively and safely remove unwanted body hair, smooth the skin, and reduce acne marks. It utilizes intelligent pulsed light irradiation technology combined with advanced optical technology for beauty care making it one of the most professional and practical products in this field.

Please read the manual carefully before use. The primary goal of the Hair Removal Device is to ensure the comfort and safety of the user. Please read all operating instructions, safety precautions, and warnings carefully before use.


IPL Principle of Operation

Using Intense Pulsed Light (IPL) technology, gentle pulses of light are applied to the skin and absorbed by the hair root. The lighter the skin and the darker the hair, the better the pulses of light are absorbed. These light pulses stimulate the hair follicle to enter a resting phase, causing the hair to shed naturally and preventing further hair growth.


Note: Use the device for 8 weeks to complete the treatment.

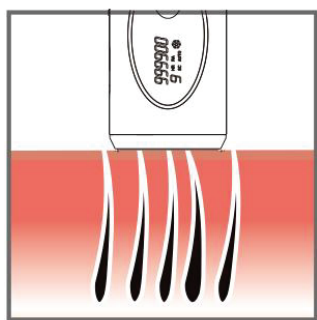

Start (1 to 4 weeks): Treat twice a week during the first 1 to 4 weeks.

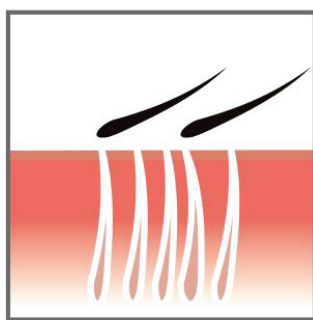

In Between (5 to 8 weeks): Continue with weekly treatments from weeks 5 to 8.

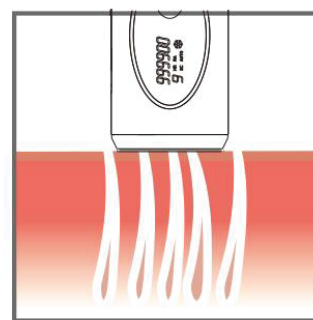

Maintenance (After 8 weeks): Continue IPL based on your needs, using the device every 1 to 2 months.


Specifications

Rated Voltage	100-240V
Rated Frequency	50/60Hz
Power Plug	UK 3-Pin Plug
Rated Power	36W
Noise Level	65dB
Light Output Area	3 cm ²
Light Source	IPL Pulse Light
Energy Level	9 Adjustable Levels
Light Output Wave Band	600 to 1200nm
Working Temperature	5 to 28°C
Light Output Time	1 to 3s
Product Net Weight	245g
Product Size	170×68×38mm

Instructions for Use




1. Shave the hair before each treatment for better results.

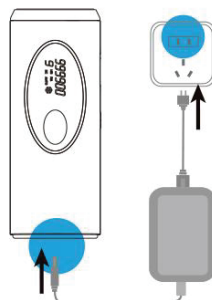


2. Connect the power supply to the device.




3. Wear goggles to protect your eyes
Note: Do not look directly at the flash.

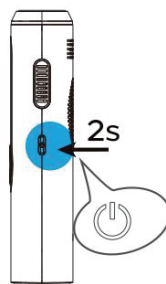


4. Long press the power button for 2 seconds to turn on the device. The energy indicator will light up and the machine will start at level 1.

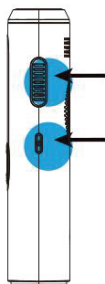


5. Short press the on/off button to adjust the energy level from 1 to 9.

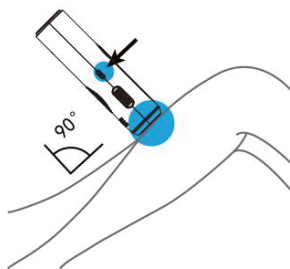

6. Hold the device vertically close to the skin, ensuring it is pressed against the skin. After the skin is detected, the display on the screen changes, and a constant flashing hair count will indicate the device is ready for use.

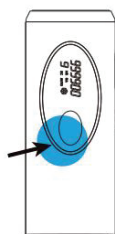


7. Press the flash button to emit the flash. To switch between auto mode and manual mode, press and hold for 3 seconds. Double-press the flash button to switch to HR/RA/SC mode.

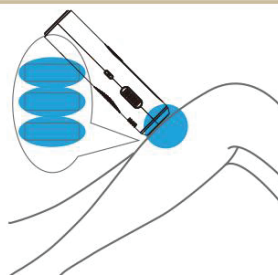

8. When using the hair removal device slide to the next area that needs to be treated. Gradually increase the level within the recommended range, and avoid flashing the same skin area multiple times.



9. Long press the "power" button for 2 seconds to turn off the device. Never unplug the device while it is still powered on to avoid potential damage.

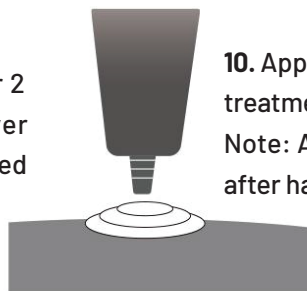

10. Apply skincare to protect your skin after treatment.

Note: Avoid sun exposure within 24 hours after hair removal.

Applicable Skin and Hair Color

Select the appropriate intensity level for your treatment. Please refer to the table to determine the intensity level that matches your skin tone and hair color, and confirm whether the product is suitable for your specific needs.


Skin Tone							
Hair Color							
White/Grey		X	X	X	X	X	X
Light Blonde		X	X	X	X	X	X
Red		X	X	X	X	X	X
Dark Blonde/ Light Brown		✓	✓	✓	✓	X	X
Brown		✓	✓	✓	✓	X	X
Dark Brown		✓	✓	✓	✓	X	X
Black		✓	✓	✓	✓	X	X

Applicable Body Parts

For Women: Ideal for head-to-toe use, including areas like the upper lip, chin underarms, legs, arms, and bikini area.

For Men: Perfect for treating the chest, back, arms, stomach, and leg hair.

Note: The recommended hair removal times mentioned above are based on research and testing. Individual hair removal times may vary from person to person.

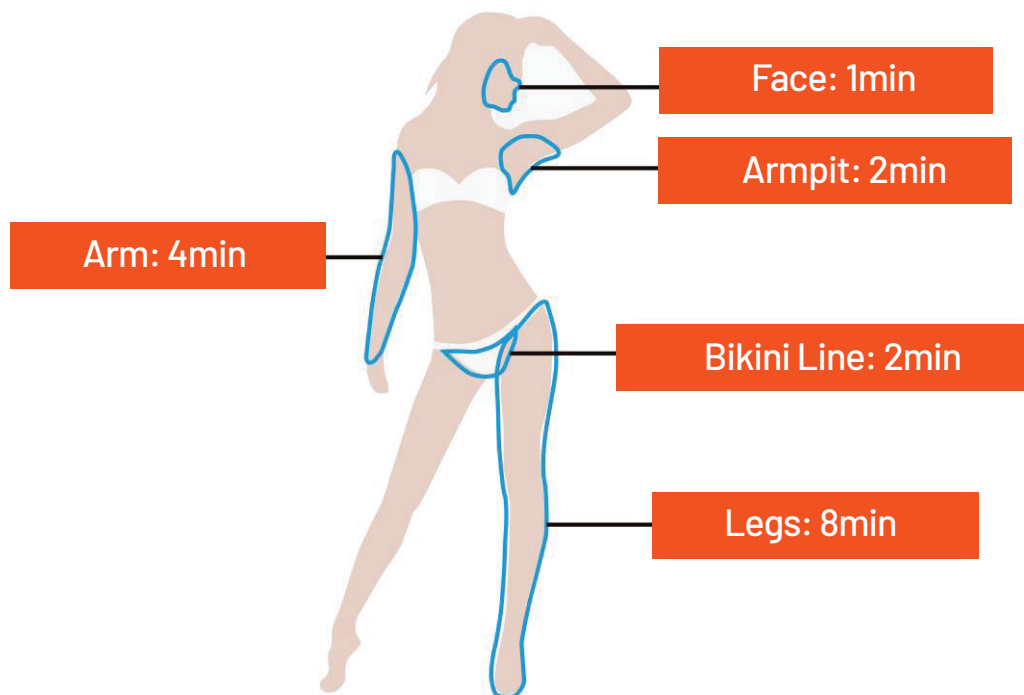

Troubleshooting

Problem	Possible Cause	Solution
The adapter or hair removal device heats up during use.	Normal Operation	No action required.
The device is turned on, but the intensity level cannot be increased or decreased.	Device needs to be reset.	To restart the hair removal device simply turn off the power and then turn it back on. If you still can't adjust the light intensity please contact the seller.
After the instrument is powered on, it can't emit light normally when placed on the skin.	The device is not positioned vertically at a 90° angle to the skin.	Once the device is positioned vertically at a 90° angle to the skin the numbers on the LED screen will flicker, indicating that the device is ready for use with normal flash mode.
After a long period of use, the numbers on the LED screen suddenly stop flickering, but the cooling fan is running normally.	The overheating protection system has been activated temporarily suspending the function. The fan continues to operate to dissipate heat.	Do not turn off the device immediately; allow it to cool down for a minute before continuing use. If it still does not flash after cooling down, try restarting it.

The numbers on the LED screen flicker, but no flash is emitted when the flash button is pressed.	Equipment system failure detected.	Restart the device. If it still does not flash after the second restart, please contact the seller for after-sales support.
Pain feeling intolerable when using this product.	You have not removed the hair before use.	If there is still hair in the area where removal is required, it may burn when using this product, causing an abnormal smell. Please ensure the hair is removed before use.
	The selected light energy level is too high.	For first-time use, it is recommended to start at the lowest gear setting and gradually increase to a comfortable level.
	This area is not suitable for using the device.	This product should not be used on private areas such as the vagina, nipples, or scrotum.
	This device is not compatible with your hair color or skin tone.	The IPL hair removal device is best suited for lighter skin tones and brown-to-dark hair. Not suitable for black skin. Not suitable for blonde, red, gray, or white hair.
The applied area became red after using the hair removal.	Slight redness is a normal reaction harmless to the human body and will fade soon.	No further action required.
No visible results after use.	The selected light intensity is too low.	Perform a skin test to determine if you can use a higher light intensity.
	For better results, follow these experimental recommendations:	Recommended usage frequency: Weeks 4–1: Use twice per week. Weeks 8–5: Use once per week. After 8 weeks: Use 2–1 times per month.
	Hair was not shaved before use.	Shave the hair on the skin's surface before each treatment for better results.

The product is not functioning properly.	This product has existing faults.	If you experience any issues, please contact the seller for assistance.
--	-----------------------------------	---

FAQ

Question	Answer
Is this product safe for use by pregnant women and children?	This hair removal device has not been tested on pregnant women or minors. Therefore, its use is not recommended for pregnant women or children.
Does this laser hair removal device work on blonde hair?	<p>How IPL Hair Removal Works:</p> <p>IPL hair removal devices emit light pulses designed to target melanin in the hair. The melanin absorbs the light generating heat that prevents hair regrowth over time.</p> <p>Without melanin, this process is ineffective. Therefore, IPL is not suitable for blonde gray, or white hair.</p>
Do I need to shave unwanted hair before every treatment, or just for the first time?	This IPL device includes a razor. For optimal results, shave unwanted hair before every treatment, not just the first time.
Is it safe to use over tattoos?	IPL hair removal should not be used over tattoos.
How often should I undergo treatment?	<p>Recommended Treatment Schedule:</p> <p>Weeks 4–1: Use twice per week.</p> <p>Weeks 8–5: Use once per week.</p> <p>After 8 weeks: Use 2–1 times per month.</p>

Is it safe to use this device for a full Brazilian treatment, including the entire bikini area?	<p>IPL Hair Removal for Bikini Area: IPL hair removal can be used for the bikini line.</p> <p>Avoid treating the genital area as it is not recommended.</p> <p>A full bikini treatment is an option, but this does not include a full Brazilian.</p>
---	--

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**