

**PORODO**  
LIFESTYLE



**Porodo Lifestyle**

**Neck and Shoulder Massager**

SKU: PD-LFST120


# Table of Contents

Precautions	2
Button Guide	3
Instructions for Use	3
Safety Warning	4
Specifications	4
Care and Maintenance	5
Troubleshooting	5
Disposal Notice	6
Warranty	7
Contact Us	7

## Precautions

- 1.** Do not pull the power cord to avoid wire breakage or poor contact.
  - 2.** When unplugging the power plug, always hold the plug and pull it out.
  - 3.** In case of any abnormal situation, such as the machine not working, stop using it immediately and turn off the power supply.
  - 4.** Do not allow infants or children to use or play with the product.
  - 5.** If you feel uncomfortable while using the product, please stop using it immediately.
  - 6.** People with the following conditions should avoid using this product: pregnancy, menstruation, high fever, malignant tumors, heart disease, recent surgery, acute illness, or other serious medical conditions.
  - 7.** This product should not be used in the following places: bathroom, kitchen, or other areas with high temperatures, high pressure, humidity, or near water and fire sources.
  - 8.** Please avoid using the product near spikes. Do not use the product if the outer cover is damaged. Contact customer service for repairs.
  - 9.** Do not use the product if the leather cover of the machine is damaged, as it may cause injury or involve the machine in a way that could lead to personal harm.
  - 10.** Always turn off the power or remove the power cord from the socket before cleaning.
  - 11.** If the power cable is damaged, please contact customer service for maintenance to prevent potential hazards.
  - 12.** It is necessary to use a 5V charger that meets national standards for charging. Otherwise, it may cause product damage or even a fire.
- Warning:** The appliance has a hot surface, and people who are sensitive to heat should take extra care when using it.

## Button Guide

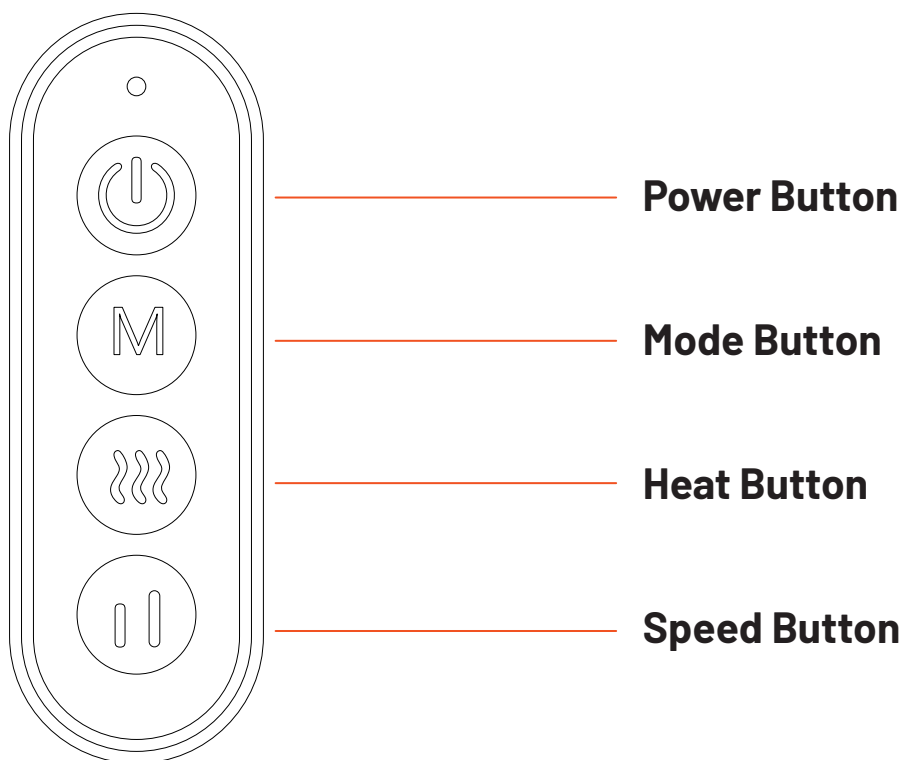

## Instructions for Use

1. The product is powered by 7.4-18650V lithium batteries. The operating time is 10 minutes, and the charging method is via 5V(2A) Type-C cable. When the product's battery is low, the red indicator will flash for 5 seconds and then automatically shut down.
2. While charging, the product cannot be used. The indicator will show a red breathing light during charging, and once charging is complete, the indicator will turn green.
3. Button Description
  - a. Power Button: Press and hold the power button for 2 seconds to start kneading and heating in the default slow mode. If the red light remains on, press and hold the button again for 2 seconds to shut down the product.
  - b. Mode Button: In the startup state, press the mode button briefly to enter the pinch mode. Press it again to return to the normal mode, and the modes will cycle in turn.

**c. Heat Button:**

**c.1** By default, the hot compress function is enabled. Tap the button to toggle between enabling and disabling the hot compress function.

**c.2** When heating is on, the red light will illuminate on the control panel. When heating is off, the green light will be illuminated.

**d. Speed Button:** Start the default slow speed by pressing the button briefly, then press quickly for a fast speed, followed by a brief press for a slow speed again. The modes will cycle in turn.

## **Safety Warning**

Stop Using the Device Immediately and Refer It to a Service Center in the Following Scenarios:

**1.** If the power cord is damaged, stop using the product immediately to avoid any potential danger.

**2.** If the fault is caused by improper use, do not remove the power plug. Please contact customer service for maintenance.

**3.** If the wire or plug of the machine is damaged, stop using the product immediately.

**4.** In the event of smoke or fire during use, immediately remove the plug from the socket and notify customer service.

**5.** If the leather or cloth parts of the product are damaged or exposed, stop using the product immediately. Contact customer service for maintenance to avoid any danger.

## **Specifications**

Rated Power	10W
Input Voltage	DC 5V/2A
Battery Capacity	2000mAh
Charging Port	Type-C

Speed Levels	2 Levels (Normal and Fast)
Charging Time	3.5 to 4 Hours
Working Time	80-120 Minutes
Working Cycle	10 Minutes
Weight	1.2kg
Product Size	178×168×1250mm
Model Number	LFS120

## Care and Maintenance

1. For daily maintenance, gently scrub the device with a neutral detergent and soft cloth. Avoid using oil, organic solvents, or other chemicals. Do not rinse the device with water.
2. When operating the product, avoid applying excessive force by pressing the switch gently.
3. Avoid heavy pressure and sharp objects that could scratch the surface of the product.
4. After use, store the device in its original box or in a dry, cool place.
5. If you find that the fault cannot be resolved during use, immediately disconnect the power to stop using the device and contact customer service. Do not attempt to disassemble the product.
6. Please unplug the power supply once charging is complete.

## Troubleshooting

**Problem:** The product does not start.

**Possible Cause:** The battery is dead.

**Solution:** Ensure the product is fully charged before turning it on.

**Problem:** The product stops while in use.

**Possible Cause 1:** The preset time has been reached.

**Solution 1:** Turn off the power before restarting.

**Possible Cause 2:** The load on the product is too heavy (e.g., stepping on or stacking heavy objects).

**Solution 2:** Reduce the load to ensure proper usage.

**Problem:** There is a sound when the product is in use.

**Possible Cause:** The sound is emitted by the mechanical structure of the product.

**Solution:** This is normal.

## Disposal Notice

The lithium-ion battery used in this product is a valuable resource that can be recycled. Please remove the lithium-ion battery and take it to a professional recycling facility for proper disposal when the product is no longer in use.



## **Warranty**

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:

**porodo.net/warranty**

## **Contact Us**

If you have any questions about this Privacy Policy, please contact us at:

**info@porodo.net**

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**