



Green Lion

Sonic Pro Rechargeable Toothbrush

SKU: GNGTB14TBBR






Model No: GTB-14

Table of Contents

| | |
|----------------------------------|---|
| Safety Precautions | 2 |
| Warning | 2 |
| Note | 3 |
| Components | 4 |
| Charging of Toothbrush | 5 |
| Charging Instructions | 5 |
| Instructions | 5 |
| Brushing Pattern | 6 |
| Function Description | 6 |
| Clean the Machine and Brush Head | 7 |
| Water Proof Function | 7 |
| Discharging the Product | 7 |
| Specifications | 8 |
| Disposal | 9 |
| Warranty | 9 |
| Contact Us | 9 |

Safety Precautions

In order to reduce the risk of failure, electric shock, injury, disaster, death, and loss of equipment or property, the following safety precautions must be observed.

| | | |
|---|---|--|
|  Danger | Potential risk of severe injury or death. |  This marking is used to remind the user of the forbidden operation steps. |
|  Warning | Potential risk of severe injury or death. |  This marking is used to remind the user of the operation steps that must be observed in order to operate the product safely. |
|  Notes | Potential risk of minor injury. | |

This product contains a built-in battery.

Do not put it into the fire for heating or charging at high temperatures, or it will cause overheating, fire, or explosion.

Warning

In case of any abnormality or malfunction, please stop using immediately and unplug the power plug.

Otherwise, electric shock or injury may be caused.

<Abnormal failure condition>

Deformation or abnormal heating of the machine or charger.

The machine or charger emits a burning smell.

The machine or charger makes abnormal sounds during use or charging.

Make sure that the working voltage of the charger is in accordance with its rated voltage, and the power plug should be inserted tightly to the root.

Otherwise, fire or electric shock may be caused.

Clean the power plug regularly to prevent dust accumulation.

Otherwise, insulation failure may be caused by moisture, thus resulting in fire.

Please pull off the power plug and wipe it with a dry cloth.

This product is not suitable for use by people with physical, sensory, or mental disabilities or lack of experience and knowledge (including children).

Guardians must supervise and guide the use of this product, and be responsible for the safety of users.

Please prevent children from treating this product as a toy.

Otherwise, accident or injury may be caused.

If the power cord of the charger is damaged, the whole charger shall be regarded as scrapped and the charger shall be replaced.

Otherwise, accident or injury may be caused.

Do not disassemble the product except when it is discarded.

Otherwise, it may cause fire, electric shock, or injury.

Do not store within the reach of children and infants.

Do not let children and infants use it.

Otherwise, accidents or injuries may occur due to careless installation of moving parts.

Do not use chargers other than those provided.

Do not use the attached charger to charge other products.

Otherwise, it may cause burns or fire due to short circuit.

Notes

If you do not use this product for a long time, please remove the power plug from the domestic power socket.

Otherwise, it may lead to electric leakage due to the decline of insulation performance, which may cause electric shock or fire.

People with severe periodontal disease, dental history, or suspected oral disease should consult a doctor before use.

Otherwise, the teeth or gums may be damaged.

Do not use this toothbrush with teeth containing a lot of friction agents, such as toothpaste for removing tobacco stains or whitening toothpaste.

People with low oral sensitivity should not use this toothbrush. Do not press the brush head against the teeth or gums, or brush a part for a long time. Do not make the machine or the plastic parts of the brush head contact the teeth.


Otherwise, it may cause tooth or gum damage.

Do not share the brush head with family members or others, as it may cause infection or inflammation.

Do not place odd objects such as coins or paper clips on the charger.

Otherwise, scale can be caused due to overheating of these chemicals. Before discarding the electric toothbrush, the battery must be removed from the toothbrush; when removing the battery, make sure it is in the power-off state; the battery should be disposed of safely.

Components

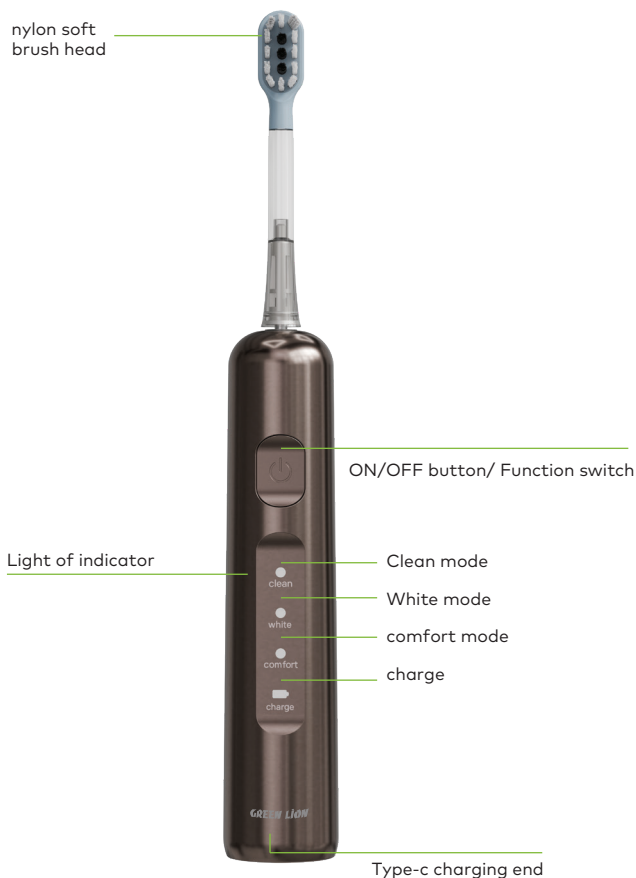

Charging of Toothbrush

If the toothbrush stops running during the initial use (just purchased or not used for 3 months or more), or if the red light of the charging indicator light flashes when using, please charge the machine.

The suitable ambient temperature for charging is 5 to 35°C. When the temperature is too high or too low, the battery cannot be charged normally or cannot be charged at all.

During use and charging, the product may become hot, but this is normal.

Charging Instructions

USB Charging: Type-C plug of charging cable is connected with the toothbrush, and USB end is connected with a 5V power supply (such as power bank, USB socket of computer, etc.).

Charging Time: 2-6 hours

Charging State: The light of indicator is flashing.

Charging Completed: The light is on when charging is full.

Low Voltage Reminder: The light flashes when using.

Instructions

1. Please insert the brush head firmly into the brush head shaft. Installing or removing the brush head in an inclined or twisted way may cause damage to the brush head or the machine.

Under the action of external force, the brush head shaft will move, which is normal during operation.

There is a small gap between the toothbrush head and the machine, which is normal.

2. Put the toothbrush in your mouth

3. Before opening, if necessary, only apply toothpaste on the brush head to prevent toothpaste splashing.


Press the gear button to select the mode

Each time the gear button is pressed, the mode is switched once.

When the brush head is pressed too hard on the teeth, the sound may increase, but this is not a fault. Please put the brush head gently against the teeth.

4. Please turn off the power supply after brushing teeth

When removing the brush head from the mouth, please do it under the condition of closed power supply, so as not to splash toothpaste foam.

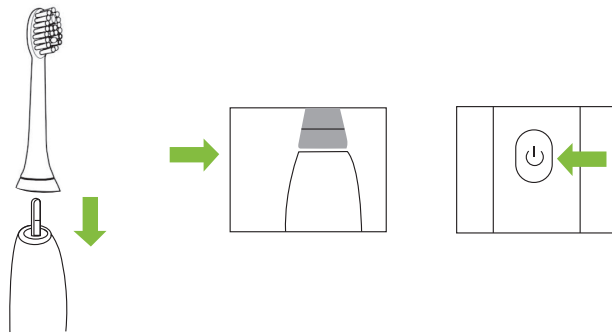

Brushing Pattern

Clean Mode: Standard mode for daily cleaning and nursing.

White Mode: Bright white mode for deep cleaning.

Comfort Mode: Low-frequency vibration mode, massages the teeth.

Function Description

Operation Mode

Sonic vibration: The highest vibration frequency is 53,400 per minute.

Timing Function


In 4 areas of the mouth, pause for 30 seconds each time and pause for 0.5 seconds in the middle.

Automatic Shutdown Function

After working for 2 minutes each time, it will shut down automatically.

Clean the Machine and Brush Head

Cleaning and maintenance must not be carried out by children. Please turn off the power supply and remove the brush head before washing with water.

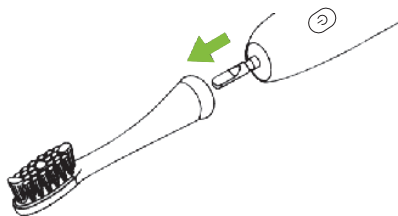


Brush Head

Rinse the brush head with running water.

Shake off the water after cleaning.

Do not remove the bristles on the brush head.

Do not wash with water higher than 80°C.

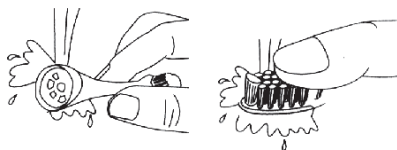

Waterproof Function

This machine has a waterproof function and cannot conduct electricity. It can be used safely in the bathroom.

Discarding of Products

The power of the product comes from recyclable lithium batteries. Please send it to the local legal recycling or disposal site, or dispose of it according to the local environmental protection regulations

Specifications

| | |
|--------------------|--------------|
| Model No. | GTB-14 |
| Material | ABS |
| Bristle Material | Nylon |
| Battery | 1200mAh |
| Power | 2W |
| Voltage | 5V |
| Vibration | 53400 VPM |
| Oscillation | 50° |
| Waterproof | IPX7 |
| Charging Port | USB-C |
| Charging Time | 2-6 hrs |
| Working Time | 120 Days |
| Brush Heads | 3 pcs |
| Product Weight | 120g |
| Product Dimensions | 210*32*29 mm |

Disposal

Do not dispose of discarded electrical appliances as unclassified urban waste. They should be handed over to an officially designated recycling center, which is beneficial for environmental protection.

Warranty

Products that you buy directly from our **Green Lion** website or shop come with a 24-month warranty.

When you buy **Green Lion** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at <https://www.greenlion.net/warranty> and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
<https://www.greenlion.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

Website: <https://www.greenlion.net/>

Sales: sales@greenlion.net

Service Support: Support@greenlion.net

Telegram: [@greenlion_global](https://t.me/greenlion_global)

Instagram: [@greenlion_global](https://www.instagram.com/greenlion_global)