



Green Lion

Sonic Plus Rechargeable Toothbrush

SKU: GNGTB15TBBK

Model No: GTB-15

Table of Contents

Specifications	2
Overview	3
BrushHead Installation	3
Safety Precautions	4
Components	5
Charging Procedure	5
Battery Indicator	5
Notes	6
Disposal	7
Warranty	7
Contact Us	7


Specifications

Model No.	GTB-15
Material	ABS
Bristle Material	Nylon
Battery	500mAh
Power	2W
Voltage	5V
Vibration	33000 VPM
Waterproof	IPX7
Charging Port	USB-DC
Charging Time	2-6 Hrs
Working Time	35 Days
Brush Heads	3 Pcs
Product Weight	57 g
Product Dimensions	208*24*24 mm

Overview



Brushhead Installation

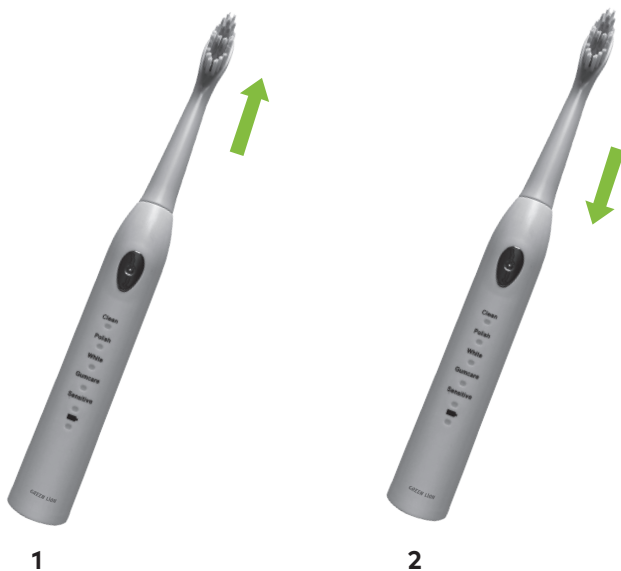

1. Removing Brush Head:

Move the brush head in the direction of the arrow.

2. Installing Brush Head:

Insert the bottom of the brush head toward the toothbrush button in the direction of the arrow.

Daily Usage Guidelines

Attach the brush head.

Wet the bristles and apply toothpaste.

Position the bristles at a 45° angle to your teeth.

Turn the brush on when it is in your mouth, to avoid splashing.

Gently glide the brush along the gum line using only light pressure from your fingertips.

Do not press or scrub hard, slowly move the brush from one area of your mouth to another.

Safety Precautions

Do not disassemble, repair, or alter the product without authorization. People who have severe periodontal disease or oral surgery should consult a doctor before use.

Not for use by or for infants and children.

Do not force the bristles hard against the teeth or gums.

Do not place the electric toothbrush in high-temperature conditions or keep it out in the sun.

Do not soak the product in water for a long time.

Please use the original charging cord provided in the box. Usage of any other accessory might cause damage to the product.

Components

Electrical Toothbrush
Brush Heads
USDC Cable
Manual

Charging Procedure

Connect the provided USB DC cable to the device.

Insert the USB cable into any USB power adapter and connect it to a power supply.

The battery symbol will light up, indicating that the brush is now charging.

Battery Indicator Description

Charge Level	The charge light flashes lowly
ON Full Charge	The charge light is on
ON Low Charge	The charge light flashes

Notes

1. Please pay attention to the status of the battery indicator during use.
2. Please charge in time when the battery is low. If you have not used the electric toothbrush for a long time, please charge the toothbrush when you use it again.



Disposal

Do not dispose of discarded electrical appliances as unclassified urban waste. They should be handed over to an officially designated recycling center, which is beneficial for environmental protection.

Warranty

Products that you buy directly from our **Green Lion** website or shop come with a 24-month warranty.

When you buy **Green Lion** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at <https://www.greenlion.net/warranty> and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
<https://www.greenlion.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

Website: <https://www.greenlion.net/>

Sales: sales@greenlion.net

Service Support: Support@greenlion.net

Telegram: [@greenlion_global](https://t.me/greenlion_global)

Instagram: [@greenlion_global](https://www.instagram.com/greenlion_global)