

PORODO



Porodo
Cordless
Heated Foot Massager
Tapping & Kneading Therapy

SKU: PDLFSTZMF026WH

Table of Contents

Product Overview	2
Specifications	2
Product Features	3
Product Schematic View	4
Package Contents	5
Operating Instructions	5
Charging and Light Indicator	6
Prohibitions	6
Safety Precautions	6
Important Notices	7
Disposal	8
Warranty	8
Contact Us	8

Overview

The Porodo Cordless Heated Foot Massager offers a soothing relief experience for tired feet, combining tapping and kneading therapy. Designed with an ergonomic foot layout, it provides an exceptional user experience with three custom modes and a heating function to enhance relaxation. This foot massager comes equipped with a 24W rated power, a 2500mAh battery capacity, and a USB-C charging port for convenience. It operates with adjustable motor speeds, ranging from 2340RPM to 2700RPM, and ensures safety with an automatic shut-off timer after two hours of use. The massager's heat temperature reaches up to 42°C, delivering a gentle warmth to further relax the feet. Additionally, the product is lightweight at 1.95kg and compact, making it easy to store and use anywhere.

Specifications

Input	DC 5V/2A
Rated Power	24W
Battery Capacity	2500mAh
Recharge Time	4H
Working Time	2H
Heat Temperature	42°C
Charging Interface	USB-C
Motor Speed	2340RPM to 2700RPM
Product Size	31.8 x 30.4 x 7.5cm
Product Weight	1.95kg
Model Number	PDLFSTZMF026

Product Features

1. Ergonomic Foot Design:

The foot massager is designed to comfortably fit the natural shape of your feet, ensuring a relaxed and therapeutic experience. Its ergonomic shape supports the feet in an optimal position for maximum comfort and effectiveness during use.

2. Three Custom Modes:

Offering three distinct modes, this foot massager provides a tailored experience to meet individual needs. Whether you prefer a light massage for relaxation, a more vigorous massage for deep relief, or a gentle heat therapy, the device allows you to select the mode that suits you best.

3. Heating Function:

Equipped with a heating function, the massager can generate heat up to 42°C, enhancing the medical effects of the massage. The warmth helps improve blood circulation, relax muscles, and provide a soothing sensation to relieve foot fatigue.

4. Adjustable Intensity Levels:

With adjustable motor speeds ranging from 2340RPM to 2700RPM, you can control the intensity of the massage. This flexibility allows you to personalize your massage experience, whether you prefer a gentle, soothing touch or a more intense, rejuvenating massage.

5. Automatic Shut-Off Timer:

To ensure safety and energy efficiency, the foot massager includes an automatic shut-off timer that turns off the device after use. This feature is ideal for those who want to relax without worrying about the device running too long.

6. Cordless Operation:

The massager operates on a rechargeable battery, providing a hassle-free, cordless experience. Enjoy a convenient and portable massage without the need for a power outlet, making it easy to use anywhere at home, in the office, or while traveling.

7. USB-C Charging:

For modern convenience, the massager charges via a USB-C port, allowing for fast and easy recharging. With a 4-hour charging time, you can quickly get your foot massager ready for the next session of relaxation.

8. Compact and Lightweight:

The foot massager is designed to be compact and easy to store, with dimensions of 7.5 × 30.4 × 31.8 cm and a weight of just 1.95 kg. Its lightweight design makes it portable and easy to carry, allowing you to enjoy foot therapy anytime and anywhere.

9. High-Quality Construction:

Made with durable materials, this foot massager is built to last, providing a reliable and long-lasting product for consistent use. Its sleek and modern design fits seamlessly into any home or office setting.


10. User-Friendly Operation:

The simple control interface ensures that anyone can operate the foot massager with ease. With intuitive buttons for mode selection and heat control, you can quickly customize your massage experience without complicated settings.

11. Battery Capacity & Recharge Time:

The massager features a 2500mAh battery capacity, providing extended use between charges. It requires approximately 4 hours to recharge fully, ensuring it's always ready when you need it.

Product Schematic View

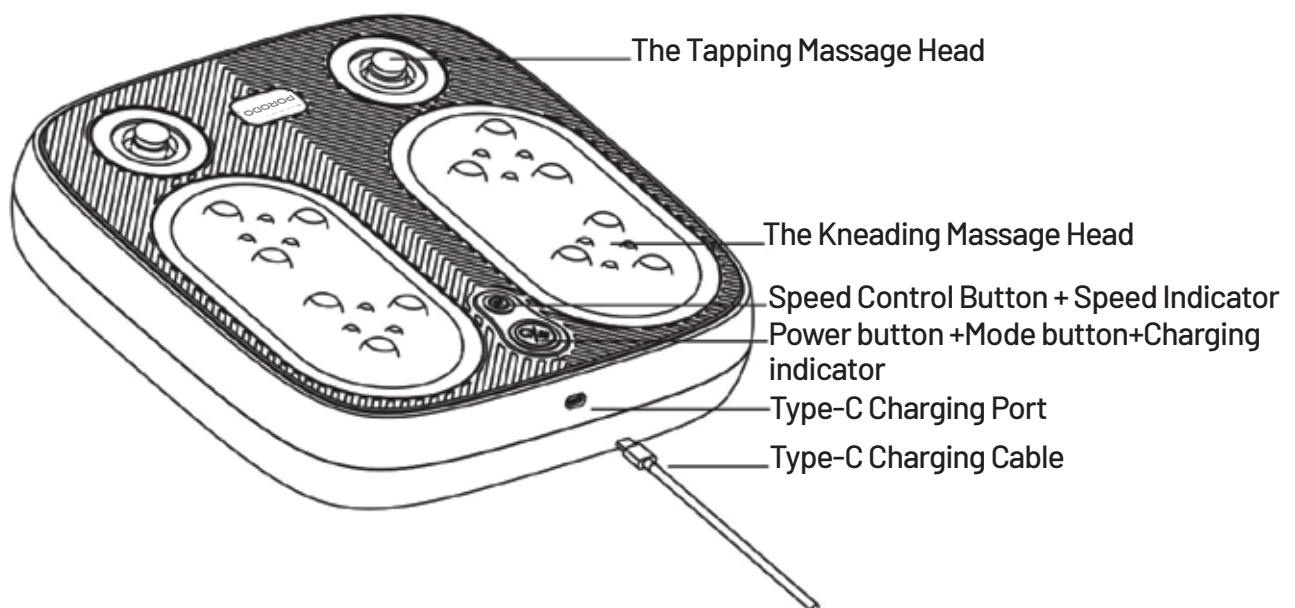

Package Contents

Host

Type-C Charging Cable

Operating Instructions

1. Power/Mode Button:

1.1. Activation: Long press to turn on the massager.

1.2. Mode Switching: Short press to cycle through modes. The blue indicator light will stay on when the device is activated, and Mode 1 is the default setting. Press once to switch to Mode 2, and press again to switch to Mode 3. This cycle repeats with each press.

1.3. Modes:

a. Mode 1: Kneading + Tapping (without heating)

b. Mode 2: Kneading + Tapping (with heating)

c. Mode 3: Kneading + Tapping (intermittent rhythm, without heating)

2. Speed Button:

2.1. Speed Adjustment: After turning the device on, press to adjust the kneading and tapping speed.

2.2. Speed Settings:

a. Speed 1 (Blue light): Default speed

b. Speed 2 (Orange light): Medium speed

c. Speed 3 (Red light): High speed

2.3. Massage Head Rotation Speed:

a. L24 :1 RPM

b. L27 :2 RPM

c. L28 :3 RPM

2.4. Tapping (Motor) Speed:

a. L2340 :1 RPM

b. L2520 :2 RPM

c. L2700 :3 RPM

3. Heating Function: To activate heating, switch to Mode 2 (red light). Heating will not work in other modes.

Charging and Light Indicator

When charging with the Type-C cable, the power mode button indicator will display a blue breathing light. Once fully charged, the indicator will show a steady blue light. If the power mode button indicator flashes, it indicates a low battery state. Please charge the product promptly.

Prohibitions

- 1.** Not suitable for cancer patients.
- 2.** Do not use directly on fracture sites for patients who have undergone metal bone needle internal fixation surgery.
- 3.** Not recommended for patients with osteoporosis.
- 4.** Should not be used by patients with venous thrombosis or pulsation.
- 5.** Not suitable for hypertensive patients.
- 6.** Pregnant women should not use this product.
- 7.** Avoid use for patients with venous dilation.
- 8.** Not recommended for patients with acute inflammation.
- 9.** Do not massage over bone areas.
- 10.** Avoid use on postoperative areas.
- 11.** Do not massage the same area for more than 20 minutes.

Safety Precautions

- 1.** Keep the massage area dry and clean. It is recommended to wear shoe covers or socks while using the product and gently press and move it.
- 2.** High-speed and high-pressure settings should only be used on the soft tissue areas of the feet, without causing pain or discomfort. Avoid using on bone areas.
- 3.** Repeated use of high-speed and high-pressure settings on the same area can

- cause abrasions, pain, or discomfort. If this occurs, stop using immediately.
4. Keep fingers, hair, or any body hair away from the massage head's moving parts to prevent pinching and injury.
 5. Do not discard or misuse this product.
 6. Do not modify the product without authorization.
 7. Do not operate or charge the product in an unattended state.
 8. This product is intended for adult use only.
 9. If you have contraindications or specific symptoms, consult a doctor before using this product.
 10. The appliance has a heating surface. Use caution if you are not sensitive to heat.
 11. Keep this product out of reach of children and pets.
 12. Avoid placing the product near heat sources, such as stoves or heaters, or in hot environments like cars.
 13. Store the product in a dry place.
 14. Do not overcharge or excessively discharge the battery, as this can damage the battery cells.
 15. Avoid using the product in environments with extreme temperatures.

Important Notices

1. Ensure there are no obstacles in front of the massage head before turning on the product. If the massage head is obstructed by external forces, the device will pause briefly and attempt to restart. Clear any blockages before resuming use.
2. This product is not waterproof. Keep it away from water or any other liquids to avoid potential damage.
3. Do not attempt to disassemble the product. If accidental disassembly occurs and causes abnormal operation, please return it for re-calibration. Any costs associated with this process will be the customer's responsibility.
4. Continuous, long-term operation may cause the motor's surface temperature to rise. Exercise caution during extended use to avoid overheating.
5. Dropping or subjecting this product to impact can cause damage and lead to malfunction. If the product is dropped or damaged, please contact customer service immediately for assistance.

Disposal

This product must not be disposed of as unsorted household waste.

It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:

info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **[porodo](https://www.instagram.com/porodo)**