

**PORODO**  
LIFESTYLE



**Porodo Lifestyle**  
**TRI-MODE HEATED**  
**Cervical Pillow Massager**

SKU: PDLFSTBZ02BK



# Table of Contents

Product Overview	2
Specifications	2
Product Design	3
Charging Instructions	4
Operating Instructions	4
Maintenance	6
Usage Guide	6
Safety Precautions	7
Disposal	8
Warranty	8
Contact Us	8

## Product Overview


The Porodo Tri-Mode Heated Cervical Pillow Massager offers a unique combination of wave, heat, and vibration features, designed to provide a relaxing and soothing experience. With three wave gears and three vibration levels, users can adjust the intensity to their preference for optimal comfort. The product is equipped with triple mode control, allowing users to customize their massage experience. Additionally, the warm wave board provides gentle heat therapy, and magnetic comfort points enhance relaxation for the neck and shoulder area. Whether you're looking for relief from tension or simply a moment of relaxation, this massager offers a versatile solution for a comforting experience.

## Specifications

Input	DC 5V
Battery Capacity	1500mAh
Working Time	1H
Recharge Time	4H
Charging Interface	USB-C
Motor Speed	3500RPM
Product Weight	590g


Material	PU + ABS
Operating Temperature	0°C~40°C
Product Size	110 x 200 x 213mm
Model Number	PDLFSTBZ02

## Product Design

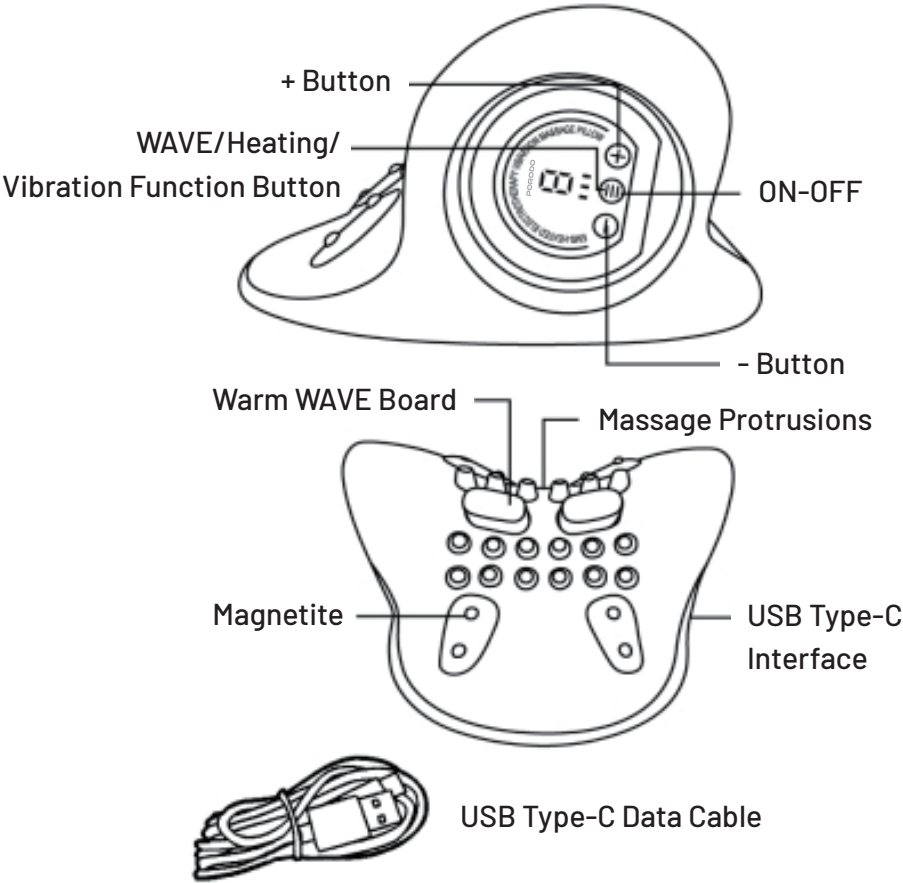


## Charging Instructions

1. Insert the included USB-C cable into the Type-C port located on the device.
2. Connect the other end of the USB cable to a USB adapter or mobile battery. Ensure the USB adapter and mobile battery are prepared separately.
3. It is recommended to use a 5V-1A power source for optimal charging.
4. While charging, three small white lights will flash. Once fully charged, the white lights will stop flashing.

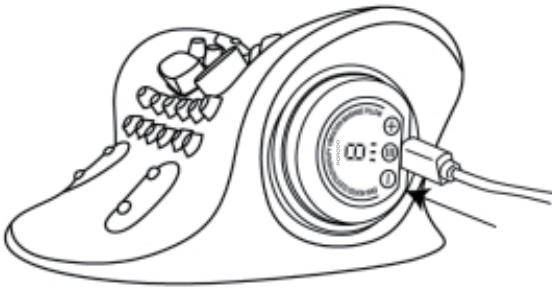

1

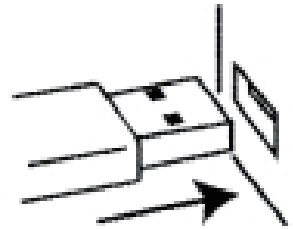

2

## Operating Instructions

### 1. Turning On the Device

Press and hold the power button for about 3 seconds to turn on the machine. The mode indicator light will appear on the first grid. To activate the WAVE function, press the increase button. When the number 1 is displayed in the gear area (out of 6), the WAVE function is enabled. Use the increase or decrease buttons to select the desired WAVE gear. When the gear area is set to 0, the EMS function will turn off.

### 2. Heating Mode

Press the power button again. When the mode indicator light appears on the second grid, the heating mode is activated. When the gear area is set to number 1, the heating function is enabled. There are 3 heating levels in total. **Warning:** Adjust heating levels based on your comfort and needs to avoid burns.

**3. Vibration Mode**


Press the power button once more. When the mode light is on the third grid, the vibration function is enabled. When the gear area is set to number 1, the vibration mode is turned on. There are 3 vibration levels in total.

**4. Automatic Shut-off**

The device will automatically turn off after 15 minutes of use. For safety, avoid using the heating function for extended periods, especially at lower temperatures (40-50°C). Prolonged use may cause skin irritation or burns. Please avoid contact with water (cooking, washing, pets, etc.) while using the WAVE function.

**5. Simultaneous Functions**

The WAVE, heating, and vibration functions can be used together or separately, depending on your preference.

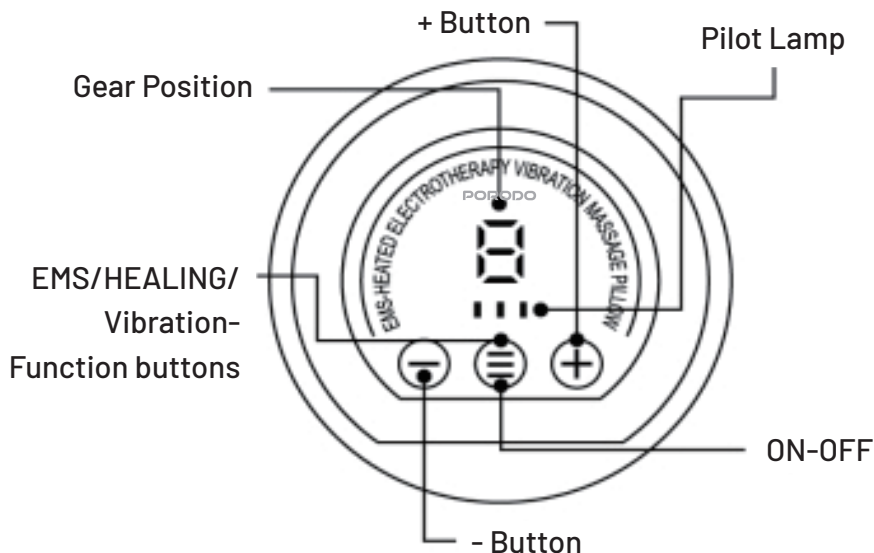

## Maintenance

1. Always use and maintain the warm WAVE board properly.
2. Wring out the cloth soaked in water, warm it up, and wipe the WAVE board thoroughly.
3. If you're concerned about dirt, soak a soft cloth with a neutral detergent diluted in warm water, squeeze it dry, and wipe off the dirt from the surface.
4. Avoid soaking the board in water or washing it all over.
5. Do not use solvents like alcohol, thinner, or benzene, as they may cause damage, discoloration, or failure of parts.

## Usage Guide

### 1. Positioning the Heater:

Adjust the position so that the hot WAVE plate is placed on the area of your neck you want to heat. The product can be used in any orientation, so feel free to try it in the direction that feels most comfortable for you.

### 2. Neck Use Only:

This product is designed specifically for use on the neck. Do not use it on other body parts, such as arms or legs.

### 3. WAVE Structure:


The WAVE plate transmits stimuli through water and moisture.

### 4. Preparation:

Before use, wet the area of your neck where the pad will touch with water or lukewarm water to activate the heating effect.

## Safety Precautions

- 1. Avoid Overcharging:** Ensure the device is fully charged before use, and do not leave it plugged in for extended periods after charging is complete.
- 2. Use Only Approved Chargers:** Only use the provided USB-C cable and charger to avoid damage to the device or battery.
- 3. Do Not Use While Charging:** For safety purposes, avoid using the device while it is charging to prevent any potential electrical hazards.
- 4. Proper Storage:** When not in use, store the device in a dry, cool place to ensure long-lasting performance and to prevent damage to the internal components.
- 5. Supervise Children:** This product is not a toy. Keep the device away from children and supervise its use to avoid any accidents or misuse.
- 6. Avoid Water Exposure:** Keep the product dry and avoid any contact with water to prevent electrical malfunctions or short circuits.
- 7. Check for Damage:** Before each use, check the product for any visible damage, including wear on the charging cable, and discontinue use if damaged.
- 8. Consult a Healthcare Professional:** If you have neck or back issues, consult with a healthcare professional before using the massager to avoid injury.

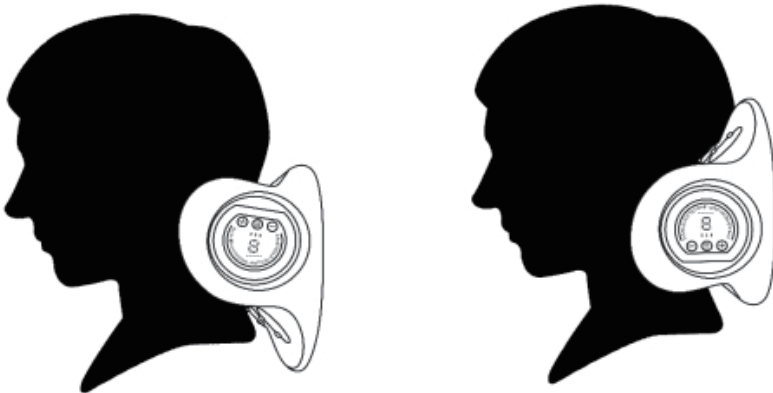



## Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

## Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:  
**porodo.net/warranty**

## Contact Us

If you have any questions about this Privacy Policy, please contact us at:  
**info@porodo.net**

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**