

# PORODO



## Porodo

### 7-IN-1 LED MICROCURRENT FACIAL MASSAGER

SKU: PDKFSTM9907WH



# Table of Contents

Product Overview	2
Specifications	2
Product Features	3
Seven Work Modes	4
Schematic View	5
Battery Status Indicator	6
Button Functions	6
Package Contents	6
Important Guidelines	6
Usage Instructions	7
Tips About Skin Care	10
Safety Precautions	10
Important Notices	11
Maintenance	11
Disposal	12
Warranty	13
Contact Us	13

## Product Overview

The Porodo 7-in-1 LED Microcurrent Facial Massager is a versatile skincare device designed to offer a range of treatments for your face. It features seven modes, including deep cleansing, red light therapy for collagen support, blue light therapy to balance oily skin, skin brightening for a radiant complexion, serum infusion to enhance product absorption, microcurrent facial lift for toning and lifting, and an eye care mode for soothing the delicate eye area. With three intensity levels and a warm massage therapy function, this device provides a customizable experience to address various skin concerns. The massager operates with a 650mAh battery, offering up to 70 minutes of use after a 3-hour recharge. It comes with a 24-month warranty, ensuring durability and quality for consistent use.

## Specifications

Input	DC 5V
Rated Power	3W
Battery Capacity	650mAh
Recharge Time	3h
Working Time	70min
Intensity Levels	3 Levels


Vibration Frequency	10000rpm
Charging Interface	USB-C
LED Light Colors	Red, Blue, Orange & Color-Mix
Product Size	50 x 54 x 165mm
Product Weight	200g
Model Number	PDKFSTM9907

## Product Features


- 1. EMS Micro-current (Galvanic Ions):** Utilizes mild electrical impulses with positive and negative ions to attract and discharge grease and grime, helping to improve skin elasticity and enhance cellular metabolism.
- 2. Vibration Technology:** Promotes the absorption of skincare products, stretches skin texture, and keeps the skin firm, supporting skin rejuvenation.
- 3. Warm Massage:** Improves the absorption of skincare products while promoting skin relaxation, enhancing texture, and maintaining skin elasticity.
- 4. Red Light Therapy (680-650nm):** Penetrates deep into the skin to improve blood circulation, reduce wrinkles, stimulate new skin cell production, and tighten the skin.
- 5. Blue Light Therapy (430-450nm):** Offers anti-inflammatory benefits, reduces sebum production, eliminates acne-causing bacteria, and improves skin appearance by enhancing oil balance.

**6. Orange Light Therapy (590-610nm):** Effectively treats red skin, relieves swelling, and promotes the metabolism of the skin's cuticles, helping to achieve a lighter, healthier complexion.


## **Seven Work Modes**

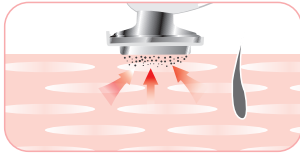


**1. Clean Mode:** The warming sensation helps open pores for deep cleaning, treating acne blemishes, and addressing dark skin spots.

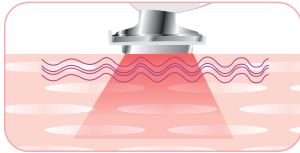


**2. RF LED Mode:** Red light activates collagen production, promoting skin elasticity, protein synthesis, and collagen formation. It works in conjunction with warm heating and vibration to smooth the skin.

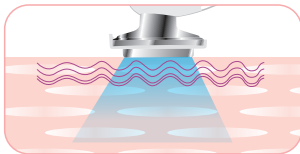


**3. Acne-T Mode:** Blue light reduces acne, effectively inhibiting inflammation and oil secretion, while improving pore size and overall skin texture.

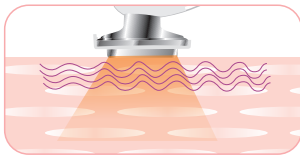


**4. Whiten Mode:** Orange light replenishes skin cell energy, accelerates metabolism, activates skin, and helps brighten the complexion.

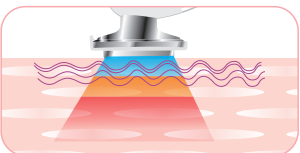

**5. Import Mode:** Microcurrent and warm heating drive skincare products deep into the skin, enhancing absorption. The combination of color light and vibration massage ensures nutrients are delivered to the skin, making it more plump and hydrated.

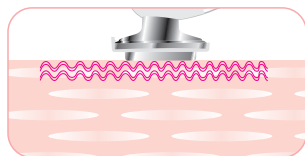


**6. EMS Up Mode:** EMS microcurrent massage stimulates blood circulation, enhances skin elasticity and firming, optimizing your skin for its best condition.

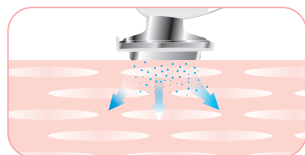

**7. Eye Care Mode:** Warm vibration massage alleviates eye fatigue, drives eye creams or essences deep into the skin, and promotes absorption, making your eyes appear refreshed and charming.

## Schematic View

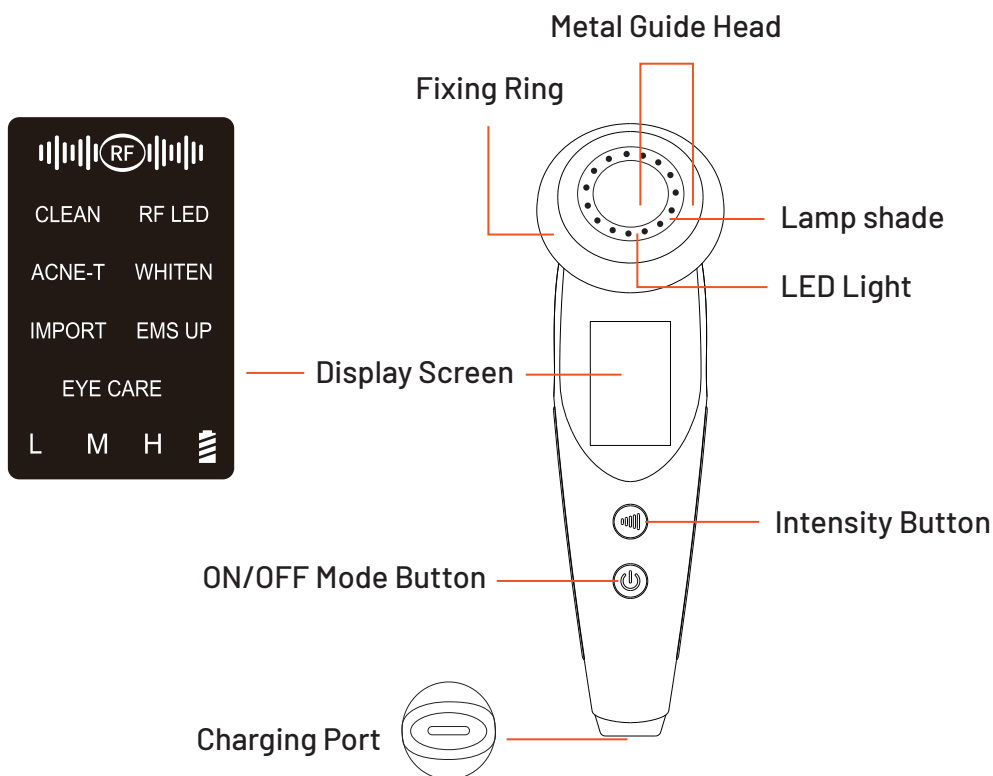

## Battery Status Indicator

White: Normal

Red: Low Battery, need charging

Blink: Charging

Bright: Fully charged

## Button Functions

### 1. Intensity ( ) button

**a.** Press briefly to change intensity.

**b.** press and hold to turn On/Off vibration in every mode.

### 2. Power ON/OFF ( ) & Mode Button

**a.** Press and hold to turn the device on and off.

**b.** Press shortly to change the mode.

## Package Contents

Main Device x1


Fixing Ring x1

Base x1

Charging Cable x1

## Important Guidelines

**1. Charge Before Use:** Ensure the product is fully charged before using it. It will automatically shut down to protect itself when charging.


2. **Skin Sensitivity:** Use the device only on the face and neck. Avoid areas with acne or redness. Sensitive skin should be used with caution. Always perform a patch test first and stop using if you experience any allergic reactions.
3. **First-Time Use:** If you are using this product for the first time or have not yet adapted to it, start with the lowest setting. If you experience any discomfort or abnormal reactions, discontinue use immediately.
4. **Skin Cleaning Step:**
  - a. Use a facial cleanser to perform the initial cleansing (without applying any product).
  - b. Select the CLEAN mode to deeply cleanse the skin. (**Note:** For the cleansing step, use makeup remover or cleansing toner, but not a cleanser or makeup-removing lotion/cleansing oil).
5. **Automatic Shut Off:** To turn the device on or off, press and hold the power button (  ). The device will automatically shut down after 5 minutes if no mode is selected.

## Usage Instructions


1. **[CLEAN] mode** (It uses the fixing ring, cotton and matches with make-up remover/cleansing toner).



a. place a cotton pad


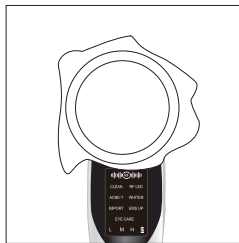


b. Fix the ring

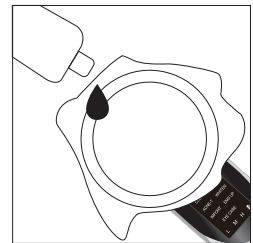

c. apply make-up remover/  
cleansing toner



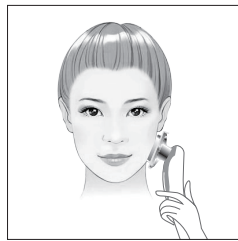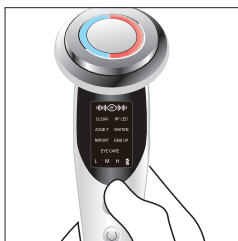


**d.** Power on: Press and hold the power button for 3s

**e.** Pick up the device

**f.** Treat your face.

**2/3/4.** [RF LED/ACNE-T/WHITEN] mode (It matches with toner / essence / lotion without fixing ring and cotton pad).


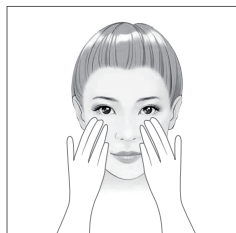


**a.** Apply toner / essence / lotion on the face evenly.

**b.** Power on: Press and hold the power button for 3s, and select your desired mode.

**c.** Hold the device, and massage your face from bottom to top.

**5.** Import mode (it matches with facial mask without fixing ring and cotton pad)

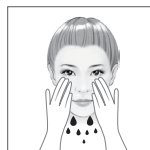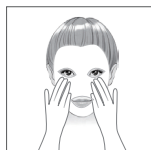


**a.** Apply facial mask on the face

**b.** Power on: Press and hold the power button for 3s and select Import mode.

**c.** Hold the device, and massage the face from bottom to top.

**d.** Remove the facial mask and clean your face.


**6. EMS UP mode** (It matches with toner/essence/lotion without fixing ring and cotton pad)

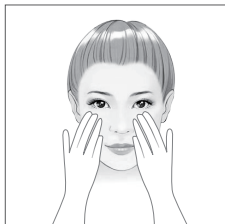

**a.** Apply lotion on the face evenly.




**b.** Power on: Press and hold the power button for 3s, and select EMS UP mode.

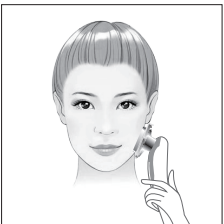

**c.** Hold the device, and massage the face from bottom to top.

**7. EYE CARE mode** (It matches with eye cream / eye essence without the fixing ring and cotton pad.

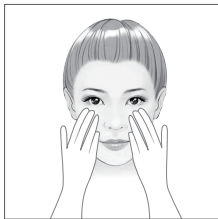


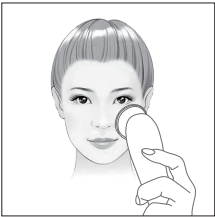
**a.** Apply eye cream or eye essence around the eye evenly.



**a.** Power on: Press and hold the power button for 3s, and select EYE CARE mode.

**c.** Hold the device, massage from inner corners to outer corners of the eye area and gently pat and press into skin.

**NOTE:** EYE CARE EMS micro-current is turned off by default, and can be turned on by pressing the (  ) button as needed.



## Tips About Skin Care

- 1.** Treatment in the morning (Preparing for make-up)
  - a.** Apply toner, essence and lotion
  - b.** Select RF LED/ACNE-T/WHITEN mode
  - c.** Make-up
- 2.** Treatment at night (Restore skin vitality)
  - a.** Face clean
  - b.** Deep clean
  - c.** Apply toner, essence and lotion (RF LED/ACNE-T/WHITEN Mode)
  - d.** Apply mask (IMPORT Mode)
  - e.** Lotion (EMS UP mode)

## Safety Precautions

- 1.** Recent facial plastic surgery or subcutaneous metal implants.
- 2.** Allergies to cosmetics or cosmetic instruments (consult your doctor before use).
- 3.** Poor health or alcohol consumption.
- 4.** Skin diseases, allergies, or fragile skin.
- 5.** Undergoing treatment for physical growth, acute illness, tuberculosis, hypertension, or infectious diseases.
- 6.** During menstruation, breastfeeding, pregnancy, or while taking sleeping medications.
- 7.** Mental health conditions or atopic skin disorders.
- 8.** Do not use on the following areas:
  - a.** Open wounds, surgical scars, oral cavity, vulva, mucous membranes, or inflamed skin areas.
  - b.** Adam's Apple, eyes, periocular area, viral warts, etc.

- c. Areas with embedded metal, silicon, or plastic implants.
- d. Recently shaved or hair removal-treated areas.

## Important Notices

1. It is recommended that the usage time does not exceed 15 minutes.
2. Do not use alcohol, thinner, or nail polish remover to clean the device.
3. Do not use the metal guide head to close the metal parts.
4. The metal guide head has a high-light surface that is prone to scratches. Please handle it with care. Accidental scratches will not affect the device's performance (this is not considered a quality issue).
5. Remove metal jewelry, glasses, contact lenses, and other metal items like necklaces, earrings, and watches before use.
6. This product is not suitable for children under 12 years old.
7. LED light can be irradiated around the eyes, so please close your eyes and avoid looking directly at the LED.
8. If you do not plan to use the product for a long time, make sure to charge it once every 6 months.
9. This product uses a lithium battery. Please dispose of it according to local regulations.

## Maintenance

### 1. Power Off the Device

Always ensure the device is powered off during maintenance to avoid any irritation or safety risks.

### 2. No Washing

Do not wash the device. Instead, wipe it gently with a dry cloth or a soft, small cloth.

### 3. Clean the Ring

If the ring becomes dirty, clean it with water. Ensure it is dry before use again.

### 4. Avoid Corrosive Solutions

Do not use alcohol, water, petrol, or nail polish remover to clean the device, as these may cause product failure, cracking, discoloration, or damage to parts.

## Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.



## **Warranty**

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **<https://www.porodo.net/warranty>** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:  
**<https://www.porodo.net/warranty>**

## **Contact Us**

If you have any questions about this Privacy Policy, please contact us at:

**[info@porodo.net](mailto:info@porodo.net)**

Website: **<https://www.porodo.net/>**

Service Support: **[support@porodo.net](mailto:support@porodo.net)**

Instagram: **porodo**