

PORODO
LIFESTYLE



Porodo Lifestyle
SMART SONIC
VISUAL ELECTRIC TOOTHBRUSH
SKU: PDLFSTETCBK



Table of Contents

Product Overview	2
Specifications	2
Product Design	3
Safety Precautions	4
App Installation	4
Connecting the App	5
Charging Instructions	5
Troubleshoot	5
Usage Guide	6
Maintenance	7
Environmental Protection	7
FCC Statement	7
Disposal	8
Warranty	9
Contact Us	9

Product Overview

The Porodo Lifestyle Smart Sonic Visual Electric Toothbrush offers an advanced oral care solution with real-time monitoring through a connected app. It features five different cleaning modes to suit a variety of needs, while the visual display allows users to track brushing progress. Equipped with IPX7 water resistance, this toothbrush is suitable for use in wet environments. It ensures efficient cleaning with advanced sonic technology and a smart brushing timer, helping users achieve optimal brushing performance. Additionally, the toothbrush provides up to 40 days of use on a full charge, making it highly convenient for daily use without frequent recharging.

Specifications

Input	DC 5V/1A
Battery Capacity	1500mAh
Water Resistance Level	IPX7
Recharge Time	3H
Charging Interface	Type-C
Working Time	40 Days on Full Charge
Vibration Frequency	40000 Stroke/Minute
Focal Length	20mm

Material	APET ABS, PP, Dupont Bristles
Camera Resolution	640 x 480 (300K)
Model Number	PDLFSTETC

Product Design



Safety Precautions

1. This product has a built-in battery. Do not put it in fire or expose it to direct sunlight, as this may damage the battery.
2. Do not use this product while charging.
3. Only use the product when the USB plug is securely in place. The product is waterproof only when the plug is properly plugged in.
4. Do not place the tooth mirror in a humid environment or liquid to avoid damage caused by immersion. Be sure to keep the plug of the tooth mirror dry.
5. When using the tooth mirror, ensure that both the socket and plug are dry before inserting the plug. Moisture may damage the host.
6. Do not attempt to modify or repair this product yourself.
7. When using the app to take pictures for diagnosis and treatment, make sure to choose a recognized dental institution.
8. Do not drop or hit the product hard.
9. Do not leave this product charging for more than 12 hours at a time. If it is not in use for an extended period, please charge it at least every six months.
10. Do not disassemble the product or replace the internal battery yourself. If you notice any swelling, leakage, or other abnormalities with the internal battery, immediately stop using the product and contact our after-sales service.

App Installation

1. Please scan the QR code below or search "SL ylz" on google play or App store and install it.



Connecting the App

1. Open the USB plug, Plug in the viewing mirror.
2. Select "Wireless LAN" Options.
3. Go to the WLAN page and find DentalCam-XXXX and connect.
4. Tap the icon to open the APP.
5. Return to the App, select "connected".
6. Enter the camera screen.

Charging Instructions

1. Press and hold the button to power off and then start charging the toothbrush.
2. A red light blinks when charging.
3. When fully charged, the green light stays on.

Troubleshoot

1. The Picture is Stuck
 - A. The data transmission distance is too far. Do not move away from the phone when using the product.
 - B. Insufficient power may also cause the product to stall. Please ensure the product is fully charged.
2. WiFi Signal Cannot Be Transmitted/Intermittently Lost
 - A. The lower the power, the weaker the WiFi signal. Ensure the product is fully charged before use.
 - B. If you are beyond the WiFi signal transmission range, try to keep the product closer to the phone.

C. Too many WiFi signals nearby may interfere with the product's signal. Try using the product in an area with fewer competing WiFi signals.


3. No Image Displayed

A. Check if your phone is connected to the correct product WiFi network.

B. After confirming the connection, if there is still no display, close both the app and the product. Then, open them again and reconnect.

Usage Guide

- 1.** Dampen the bristles of the toothbrush and apply some toothpaste to them.
- 2.** Press the On/Off button to turn on the toothbrush.
- 3.** Gently apply light pressure while brushing to maximize the effect of the toothbrush.
- 4.** Slowly move the toothbrush at a 45° angle between the bristles and gums.

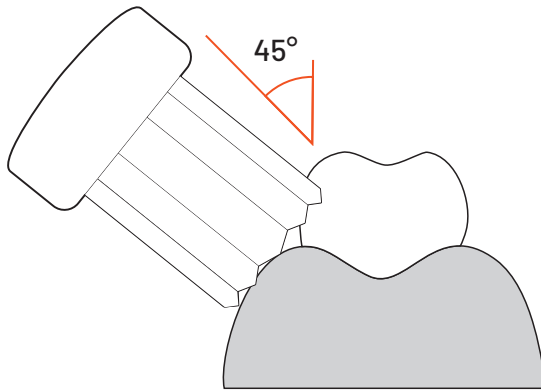

- 5.** During brushing, press the button to switch between brushing modes as needed.

NOTE: The recommended brushing time is 2 minutes.

- 6.** To power off, press and hold the button for more than 5 seconds.

- 7.** After brushing your teeth, clean the toothbrush head with water (make sure to cover the charging bottom when charging).

Maintenance

1. Wash the brush head and bristles after each use and allow them to air dry naturally. Please do not pull the bristles, as this may cause them to deform and fall off.
2. Wipe the brush handle with a dry cloth after each use and let it air dry naturally.
3. Do not use water hotter than 40°C or acidic cleaning solutions to clean the brush handle, charging cable, power cord, or adapter. The adapter should never be submerged in water.
4. Before cleaning the charging device, ensure that the adapter is unplugged from the socket to avoid the risk of electric shock.
5. It is recommended to replace the brush head every three months, or sooner if the bristles are bent or deformed.

Environmental Protection

When disposing of products, do not place them with general household waste. Built-in rechargeable batteries contain substances that may harm the environment. When discarding products, do not throw them away carelessly. Please ensure they are disposed of at a government-designated recycling station.

FCC Statement

1. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference; and

(2) This device must accept any interference received, including interference that may cause undesired operation.

2. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for assistance.

RF warning statement:

The device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure conditions without restriction.

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.

Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **<https://www.porodo.net/warranty>** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
<https://www.porodo.net/warranty>

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **<https://www.porodo.net/>**

Service Support: **support@porodo.net**

Instagram: **[porodo](#)**