

PORODO
LIFESTYLE



Porodo Lifestyle

**RED LIGHT THERAPY
SCALP MASSAGE COMB**

SKU: PDLFSTFH918BK

Table of Contents

Product Overview	2
Specifications	2
Product Design	3
Package Contents	3
Safety Precautions	3
Usage Guides	4
Liquid Addition	4
Timing Function	5
Power Off	5
Cleaning Guide	5
Usage Warnings	5
Charging Instructions	6
Maintenance and Storage	6
Troubleshooting	7
Battery Removal and Disposal	7
Disposal	8
Warranty	9
Contact Us	9


Product Overview

The Porodo Red Light Therapy Scalp Massage Comb offers a unique and effective solution for scalp care. It features customizable massage modes, allowing you to personalize your experience. The 3ml liquid tank ensures efficient oil application, while the 360° scalp massage enhances relaxation and promotes blood circulation. With its built-in red light therapy, this device helps improve scalp health. It also includes a safety auto shut-off feature for added peace of mind. This comb is a perfect addition to your self-care routine, offering a combination of comfort and therapeutic benefits.

Specifications

Input	DC 5V
Rated Power	10W
Liquid Tank Capacity	3mL
Battery Capacity	1400mAh
Charging Interface	Type-C
Recharge Time	3H
Product Size	50 x 85 x 212mm
Product Weight	250g
Working Time	40 Minutes
Model Number	PDLFSTFH918

Product Design

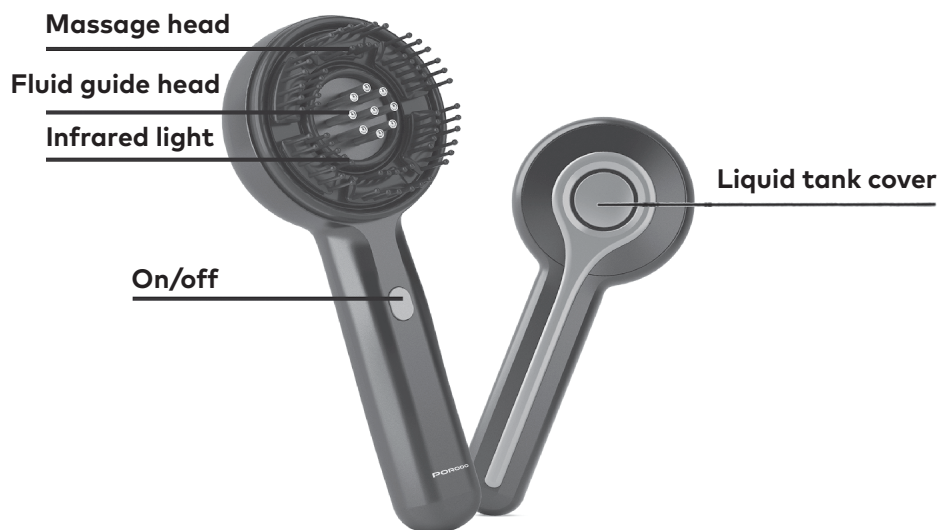

Package Contents



Safety Precautions

1. This product is intended for household use only and is not a substitute for medical devices. It is not suitable for diagnosis, treatment, or commercial purposes.
2. Before use, remove any accessories or items that may interfere with the massage function.
3. Discontinue use immediately if you experience any pain or discomfort.

4. If an abnormal issue arises during use, turn off the power and contact customer service for assistance. Do not attempt to disassemble or repair the product yourself.
5. When using in specific environments, always comply with relevant safety regulations.
6. Use only charging equipment that matches the product's input voltage. Ensure safe electricity consumption and avoid mixing with other charging devices.
7. Keep the product away from fire or corrosive environments.

Usage Guides

1. Charging:

Do not use the product while charging. Once charging is complete, unplug the power adapter and proceed to use the product.

2. Boot/Mode Selection:

- a. Press and hold the on/off button for about 2 seconds. A prompt tone ("beep") will indicate the product has powered on. The default mode is Mode 1.
- b. After the device powers on, press the on/off button to switch between modes.
- c. The prompt tone "Drip" confirms a successful switch.

3. Modes:

- a. Mode 1: Quick massage + infrared light
- b. Mode 2: Slow massage + infrared light
- c. Mode 3: Slow and fast massage + infrared light

Liquid Addition

1. Open the liquid tank cover and add the appropriate product according to your desired dosage. It is recommended to add less than 3ml to avoid wastage.
2. Clean the liquid storage bin after each use to maintain product cleanliness.
3. Different brands and concentrations of care solutions may affect the ball discharge speed. A lower concentration may cause minor leakage, so be sure

to clean up promptly.

4. Special attention: Do not add unknown substances or corrosive liquids to avoid damage to the product.

Timing Function

Smart timing for 10 minutes

Power Off

Long press the switch button to shut down the device.

Cleaning Guide

1. Appearance Cleaning:

If needed, rinse the plastic surface with water directly. Avoid getting water in the bottom charging port.

2. Liquid Storage and Outlet Cleaning:

After each use, rinse the liquid storage bin and outlet hole with clean water to prevent blockage. Ensure the product is wiped clean and placed in a ventilated, dry location.

Usage Warnings

1. During normal operation, the product may emit sounds such as "hum" or "click." These are typical and indicate the product is functioning as intended.

2. If the product is new and unsealed, there may be a small amount of liquid in the storage bin. This is a result of the factory's quality inspection process, and it is safe to use.

Charging Instructions

1. Connect the Type-C charging cable to the charging port of the device.
2. Plug the other end of the charging cable into a USB port or 5V adapter and connect it to a power socket to charge the device.
3. When the device is low on power, the indicator will flash. Once charging is complete, the red light will remain on.
4. The charging time is approximately 3 hours.

Important Notes:

- a. Ensure the charging port is dry before charging.
- b. It is normal for the device and power adapter to heat up slightly during charging.
- c. At the End of Charging, unplug the power adapter and disconnect the charging cable from the device.

Maintenance and Storage

1. Storage

- a. If the storage period exceeds 3 months, recharge the machine before use.
- b. Avoid placing the machine on the floor, walls, or roof.
- c. Ensure the storage room is well-ventilated to prevent the buildup of strong sunlight and corrosive gases.

2. Transportation

- a. Use product packaging with shock-proof features suitable for air, railway, highway, and sea transport.
- b. Protect the machine from rain, snow, splash, inversion, and collision during transportation.

3. Power Check

Verify the power supply. If the voltage is outside the specified range, do not turn on the machine.

4. Upkeep

- a. Prevent moisture, rodent activity, and pollution by taking appropriate

measures. Clean the machine while it is powered off.

b. After use, place the instrument in a safe location and avoid placing heavy objects on it.

5. Failure and Repairs

In case of failure, contact professional personnel for repairs. Non-professional personnel should not attempt to disassemble or assemble the machine.

Troubleshooting

If a fault occurs during operation, please follow the steps below. If the issue persists or the product is still not functioning properly, disconnect the power supply. The product should only be repaired by qualified technicians. Please avoid disassembling the product.

Problem	Possible Cause	Solution
Unable to boot	Check to see if the battery level is low	If the battery level is low, charge it
Weak massage	Check if the battery is weak	If the battery level is low, charge it


Battery Removal and Disposal

WARNING: Non-professionals should not attempt to disassemble the product themselves.

1. The product must be powered off before removing the battery.
2. Disconnect the external power supply and remove the product casing.
3. Unplug the battery.
4. Remove the battery.
5. Place the used battery in a designated recycling location.
6. Follow local environmental protection laws and disposal guidelines.

NOTE: Before Disposal, ensure the battery is removed according to local regulations. Dispose of the battery in compliance with applicable laws and regulations.

Cover Opening Guide:

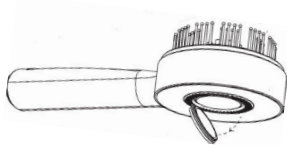

Disposal

This product must not be disposed of as unsorted household waste. It is important to separate such waste for proper treatment and recycling, in compliance with local waste management regulations.



Warranty

Products that you buy directly from our **Porodo** website or shop come with a 24-month warranty.

When you buy **Porodo** products from any of our approved sellers, you only get a 12-month warranty. If you want to extend this warranty, go to our website at **porodo.net/warranty** and fill out the form with your information. Don't forget to upload a picture of the product too. After we've checked and accepted your request, we'll send you an email to confirm that your product's warranty has been extended.

For more info, please check:
porodo.net/warranty

Contact Us

If you have any questions about this Privacy Policy, please contact us at:
info@porodo.net

Website: **porodo.net**

Service Support: **support@porodo.net**

Instagram: **porodo**